

(19)



Europäisches Patentamt
European Patent Office
Office européen des brevets

(11) Publication number:

0 327 209
A2

(12)

EUROPEAN PATENT APPLICATION

(21) Application number: 89300296.4

(51) Int. Cl. 4: **A61M 5/14**

(22) Date of filing: 13.01.89


(30) Priority: 15.01.88 US 144795

(43) Date of publication of application:
09.08.89 Bulletin 89/32(84) Designated Contracting States:
AT BE CH DE ES FR GB GR IT LI LU NL SE(71) Applicant: **SHERWOOD MEDICAL COMPANY**
1831 Olive Street
St. Louis, MO 63103(US)(72) Inventor: **Jackson, Edward**
18 Gail Drive
Northford Connecticut 06472(US)
Inventor: **Pasqualucci, Joseph**
12 Spruce Brook Road
Seymour Connecticut 06483(US)(74) Representative: **Brown, Keith John Symons et al**
c/o John Wyeth & Brother Limited Patent
and Trademark Department Huntercombe
Lane South
Taplow Maidenhead Berkshire SL6 0PH.(GB)

(54) Motor unit for a fluid pump and method of operation.

(57) An enteral nutrition pump system (10) operates in a cyclical manner with a period between cycles being selected in accordance with the desired fluid delivery rate. Each pump cycle may correspond to a single rotation of the rotor (18) or a fractional rotation of the rotor. Rotor rotation may alternatively be sensed by utilization of magnetic sensors (50) or by monitoring of the AC component of current supplied to a DC motor driving the rotor.

FIG.1

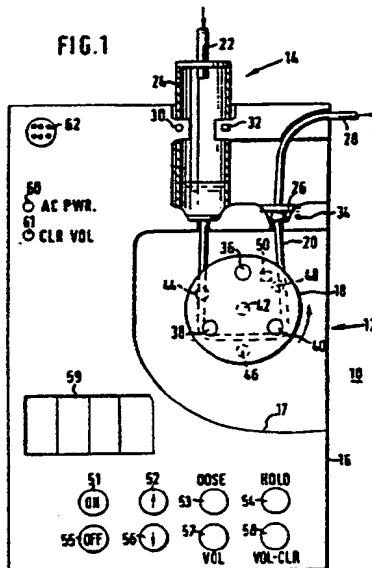

Xerox Copy Centre

EP 0 327 209 A2

BEST AVAILABLE COPY

MOTOR UNIT FOR A FLUID PUMP AND METHOD OF OPERATION

This invention relates to pumps for delivering medical fluids and particularly relates to peristaltic pumps for delivery of enteral nutrition fluids to a patient.

In accordance with known techniques the delivery of enteral nutrition fluids to a patient can be accurately controlled as to volumetric delivery rate by the use of a delivery system which includes a motor unit and a disposable delivery set. Likewise similar systems may be used for pumping of other fluids for medical purposes, such as intravenous infusion, blood pumping or supply of measured volumes of liquid medication to pre-loaded syringes or other containers.

In known systems for delivering enteral fluids the rate of fluid delivery is controlled by regulating the speed of a pump motor in accordance with the desired volume rate. Pump motor speed may be controlled, for example, by providing pulses to a stepper motor. Another system for providing variable rate fluid delivery makes use of a peristaltic pump with variable tension on the pump tube in combination with a constant speed motor.

In other known systems for pumping medical fluids there are provided means for monitoring rotation of the pump rotor, for example, by magnetic detection or by optical rotation detectors. In such systems the actual rotation rate of the motor is compared to the desired rotation rate for purposes of making corrections to the rotation rate of the motor. Alternately the motor may be operated to rotate the pump by a number of rotations corresponding to the desired volume.

It is an object of the present invention to provide a new and improved method for regulating the volumetric rate of fluid delivery in a medical fluid delivery system and to provide apparatus for carrying out the improved method.

In accordance with the present invention there is provided a motor unit for a medical fluid delivery system for use with a disposable delivery set for pumping medical fluid, characterised in that it comprises: pump operating means, including a motor, for acting in cooperation with the delivery set to deliver a volume of the fluid during each operating cycle;

and pump control means for controlling the pump operating means to deliver the fluid at a desired volumetric rate, the pump control means including means responsive to signals indicating the desired rate for activating the pump operating means for one of the operating cycles and for repeating the activation at variable time intervals selected in accordance with the desired volumetric rate.

In accordance with a preferred embodiment of the invention the pump operating means is a pump rotor for operating in connection with a pump tube on the delivery set to form a peristaltic pump and the pump operating cycle comprises a selected angular rotation of the pump rotor. In one arrangement the pump control means includes means for sensing the condition of the pump operating means with respect to an operating cycle, and the sensing means comprise a magnet and a magnetic field sensor. In another arrangement the sensing means may detect the AC component of the current supplied to a DC motor.

The medical fluid delivery system which comprises the novel motor unit and a disposable fluid delivery set carries out a novel method for controlling the rate of fluid delivery. The novel method includes providing means for detecting the completion of an operating cycle of the pump operating means and operating the motor unit until the completion of the operating cycle is detected. Operation of the motor unit is repeated at variable time intervals which are selected in accordance with the desired rate of fluid delivery.

In accordance with another aspect of the present invention the AC component of the DC motor current is detected and compared to a reference level in order to detect the current variation which results from the presence of a pump tube on the rotor. Accordingly, any mis-installation of the pump tube will be detected by the fluid delivery pump.

For a better understanding of the present invention, together with other and further objects, reference is made to the following description, taken in conjunction with the accompanying drawings, in which:

Figure 1 is a plan elevation view of an enteral fluid delivery system incorporating the present invention.

Figure 2 (comprising Figures 2a, 2b and 2c) is a circuit diagram for the system of Figure 1.

Figure 3 is a circuit diagram for a portion of a modified delivery system in accordance with the present invention.

Figure 4 is a timing diagram illustrating signals utilized in the present invention.

Figure 1 is an illustration of an enteral fluid delivery system incorporating a motor unit in accordance with the present invention. The enteral delivery system 10 includes a motor unit 12 and a disposable delivery set generally indicated as 14 which is arranged to be mounted on the motor unit. The motor unit 12

includes a housing 16, which in the illustrated embodiment includes a recess 17 within which a rotor 18 is mounted. Rotor 18 is driven by a conventional constant speed DC motor which drives shaft 42. The delivery set 14 includes a pump tube 20, made of flexible plastic which surrounds rotor 18 and interacts with 3 rollers 36, 38 and 40 mounted on rotor 18 to form a peristaltic pump. Rotation of the rotor 18 in the direction indicated by the arrow in Figure 1 causes the rollers 36, 38 and 40 to interact with pump tube 20 and pump fluid through the tube at a rate which is determined by the rate rotation of rotor 18.

Delivery set 14 includes an inlet tube 22, which is connected to a supply of enteral fluids, such as a fluid reservoir which may be mounted on an IV pole above the motor unit 12. The inlet tube 22 is connected to drip chamber 24 which is mounted to a recess on housing 16 and secured to one end of pump tube 20. The outlet end of pump tube 20 is provided with a mounting member 26 which is received in another recess on housing 16 to thereby secure the outlet end of tube 20. A fluid delivery tube 28 is connected to mounting member 26 and supplies fluid pumped by the system to an enteral feeding tube connected to a patient or to another medical fluid delivery system.

The system illustrated in Figure 1 additionally includes a light source 30 and a light detector 32 for operation in connection with drip chamber 24 to detect the occurrence of drops in the drip chamber which pass between light source 30 and detector 32 in a manner which is known in the art. Mounting member 26 includes magnetized material, the presence of which can be detected by magnetic field detector 34.

The motor unit 12 includes control buttons 51 through 58 for operating the unit to turn it on or off, to set the dose or volume rate of fluid delivery by the pump, to interrupt operation of the pump and to increase or decrease the designated fluid volume or volume rate. A four digit alphanumeric segment display 59 is provided for indicating the selected fluid delivery rate or delivered volume and for providing alarm messages or codes. Light emitting diode 60 and 61 are provided for indicating that the unit is plugged into AC power or indicating that the volume setting has been cleared. An enunciator 62 is provided for signalling an audible alarm to indicate, for example, that the pump has completed delivering a designated volume of fluid.

Housing 16 is provided with a magnetic field sensor 50 which is arranged adjacent and behind rotor 18 in order to detect the magnetic field provided by magnets 44, 46 and 48 which are mounted on rotor 18. The presence of the magnets 44, 46 and 48 is detected as the magnets pass sensor 50 during rotation of rotor 18.

Figure 2 is a schematic diagram of the circuits in the pump motor unit 12 of Figure 1. The schematic representations of the various components of Figure 1 have been given same reference numerals in Figure 2.

The motor unit operates under the control of a microcomputer 64 which is provided with a control program which is set forth in Appendix I. A programmable interval timer 68 is provided for operating and initiating microcomputer 64. A clock 66, operating at 2 Mhz, provides clock pulses to the system. The various controls of the unit, 51 through 58, are provided as input signals which ground various input terminals of the microcomputer to thereby signal the operators input instructions. The alphanumeric display 59 is driven by the microcomputer as is LED indicator 61. Additional inputs to the microcomputer are provided by the magnetic field sensors 34 and 50 which sense respectively the magnetized mounting member 26 and the magnets 44, 46 and 48 on rotor 18. Likewise the drop detector 30, 32 is connected to provide input signals to the microcomputer. An AC power rectifier 72 is provided for AC operation and battery charging. Portable DC operation is available using battery 74. The AC circuit is arranged to charge the DC battery when the unit is connected to AC power. A low battery and dead battery detector circuit 70 is provided to signal the microcomputer that the battery needs recharging. The microcomputer provides an output motor signal which is coupled by transistor 80 to switching transistors 82, 84. Transistor 82 turns on the power supply to motor voltage regulator 88 when the motor is to be operated and transistor 84 short circuits the motor to lock it into position when the motor signal is no longer present. Switching transistor 78, which is provided with a power signal by transistor 76, operates to supply current to the motor system and the other electronic systems by voltage regulator 86 when power is turned on. The motor 90 is provided with a safety circuit 92 which provides a short circuit when the motor is operated for an excess period of time. The short circuit causes fuse 94 to open thereby disabling the set when continuous motor operation occurs, to avoid providing excess enteral fluid to a patient.

Unlike conventional enteral nutrition systems the system 10 of the present invention is designed to provide an intermittent motor operation with the periodicity of the intermittent operation being regulated to adjust to the desired rate of fluid delivery. The operation of the system of the present invention is therefore cyclical and will be explained with respect to timing diagrams of Figure 4. Graph A of Figure 4 illustrates the motor voltage of the enteral fluid delivery system 10. The motor voltage is turned on and operated for a time period G which is regulated by detecting the rotation of rotor 18, in the case of Graph A for one

complete revolution. With reference to Figure 1 it may be seen that during one complete revolution, represented by motor voltage period G, three magnets 44, 46 and 48 all pass magnetic field detector 50 and are sensed thereby. Curve B in Figure 4 illustrates the output signal from the rotor sensing magnetic field detector 50 which occurs during the cycle of operation indicated by motor voltage G. During an initial period of approximately .45 seconds designated F in Figure 4 the operation of the rotor sensing is inhibiting by software in microcomputer 64, so that the initial on period J of magnetic field detector 50 is not responded to by the program. Thereafter, during one complete revolution of the rotor, the signal from detector 50 goes to zero as each magnet is encountered by detector 50. Upon detection of the third magnet, at the end of period G, the motor voltage is turned off. In accordance with the preferred embodiment of the present invention the unit repeats the cyclical operation a time period I after initiation of the first operation. The time period H during which there is provided no motor voltage is permitted to be variable, since it depends on the actual time taken for rotation of the rotor and the selected interval I. The interval I is selected according to the rate of fluid delivery to be provided by the set which is set by the operator. In one embodiment of the invention period I varies from 13.5 seconds corresponding to a delivery rate of 100 milliliters per hour to 4.5 seconds corresponding to a fluid delivery rate of 300 milliliters per hour. Motor operation period G takes approximately 4 seconds but may vary according to mechanical conditions of the motor and pump tube.

Graph E in Figure 4 shows an alternate timing arrangement wherein the motor cycle consists of a single one-third of a rotation of the rotor 18. In accordance with the operation method of Graph E the motor current period G' is ended by the detection of the first of the three circumferentially arranged magnets by magnetic field sensor 50. Again the timing I between each operating cycle of the motor is varied in order to control volumetric fluid rate delivered by the pump. In the same embodiment as previously discussed a fluid rate of 1 to 100 milliliters per hour can be delivered using a cycle interval I which ranges from 450 to 4.50 seconds.

As an alternate, or in addition to providing magnets on rotor 18 for purposes of detecting completion of a motor cycle, the motor current may be monitored for purposes of determining the rotational position of the rotor 18. Figure 3 is a schematic diagram of a circuit wherein there is provided a motor current monitoring circuit 96 which includes a low resistance resistor in series with the motor the voltage across which is AC coupled to an AC amplifier 92 for purposes of monitoring the AC component of the DC motor current. Graph C of Figure 4 illustrates a typical motor current for the operating cycle of Graph A of Figure 4. The motor current initially rises to a high level for purposes of overcoming the starting resistance and accelerating rotor 18 to its normal velocity. Thereafter the motor current drops but reaches periodic peaks corresponding to the resistance of rollers 36, 38 and 40 as they stretch pump tube 20 to its furthest position. The peak periods of motor current, which are illustrated as negative going pulses in the digitized signal of curve D, which is the output at point 96 of the circuit of Figure 3, may be used for purposes of detecting rotor position and may be used also for assuring that the pump tube 20 is properly mounted to rotor 18. Because of the initially high rotor current, which results from starting up the rotor, the current sensing is software inhibited for time period K of approximately .25 seconds prior to initiating the threshold detection which results in the pulses of curve D. Each of the pulses illustrated in curve D, which are negative going, have a positive going pulse which occurs a time period L prior to the end of a motor cycle, there being three such pulses during one rotation of the rotor. Accordingly, the curve D signal can be used for purposes of detecting and monitoring rotation position of rotor 18, and thereby indicating to the motor control circuit the completion of an operating cycle. As an alternate to providing delay L after the end of the curve D pulses, the motor cycle may be arranged to end at the end of the pulse, providing a different rotor position between cycles.

The motor current monitoring previously described can additionally be used in cases wherein the motor voltage is provided only for a one-third rotation of the rotor as discussed with respect to curve E.

An additional use of the motor current monitoring circuit, which provides the signal of curve D of Figure 4 is to provide assurance to the system that the pump tube 20 has been properly mounted on rotor 18. Accordingly at the initiation of motor current and after a delay period K a flag can be set by the microprocessor which is cleared by the negative going pulse of curve D to indicate proper pump tube positioning. The flag would be reset at the start of each operating cycle or may also be reset on the occurrence of the one-third rotation of the rotor sensing current shown by curve B. If the flag is not cleared by the negative going pulse of curve D there is an indication that either there is no pump tube or that the pump tube has been improperly mounted and an alarm signal can be initiated.

Claims

1. A motor unit for a medical fluid delivery system for use with a disposable delivery set for pumping medical fluid, characterised in that it comprises:
 - 5 pump operating means, including a motor, for acting in cooperation with the delivery set (14) to deliver a volume of the fluid during each operating cycle;
 - and pump control means for controlling the pump operating means to deliver the fluid at a desired volumetric rate, the pump control means including means responsive to signals indicating the desired rate for activating the pump operating means for one of the operating cycles and for repeating the activation at
 - 10 variable time intervals selected in accordance with the desired volumetric rate.
2. A motor unit as claimed in claim 1 wherein the pump operating means comprises a pump motor (18) for acting in cooperation with a pump tube (20) on the delivery set (14) and wherein the operating cycle comprises a selected angular rotation of the pump rotor.
3. A motor unit as claimed in claim 1 or 2 wherein the pump control means includes means for sensing
 - 15 the condition of the pump operating means with respect to the operating cycle.
4. A motor unit as claimed in claim 3 wherein the means for sensing comprises at least one magnet (44, 46, 48) and a magnetic field sensor (50).
5. A motor unit as claimed in claim 3 wherein the pump operating means comprises a DC motor and wherein the means for sensing the condition of the pump operating means comprises means for sending
 - 20 the AC component of the current supplied to the DC motor.
6. A motor unit for a medical fluid delivery system for use with a disposable delivery set having a flexible pump tube for engagement with the motor unit for pumping medical fluids, characterised in that it comprises:
 - a housing (16) for receiving at least portions of the delivery set (14);
 - 25 a motor driven rotor (18) mounted on the housing (16) for rotational movement with respect thereto and for receiving the flexible pump tube (20) for forming a peristaltic pump;
 - at least one magnet (44, 46, 48) and at least one magnetic field detector (50) mounted to the rotor (18) and the housing (16) for relative rotational motion with respect to each other, the magnet and the magnetic field detector being operatively adjacent to each other at at least one rotational position of the rotor;
 - 30 and pump control means, responsive to selection of a fluid delivery rate and responsive to the magnetic field detector (50), for operating the motor until the rotor is in the rotation position with the magnet and the magnetic field detector adjacent each other and for repeating the operation at variable time intervals selected in accordance with the selected fluid delivery rate.
7. A motor unit for a medical fluid delivery system for use with a disposable delivery set having a
 - 35 flexible pump tube for engagement with the motor unit for pumping medical fluids, characterised in that it comprises:
 - a housing for receiving at least portions of the delivery set;
 - a rotor driven by a DC motor, mounted on the housing for rotational motion with respect thereto and for receiving the flexible pump tube for forming a peristaltic pump;
 - 40 means for sensing the AC component of the current supplied to the DC motor thereby to sense rotation of the rotor by an incremental amount corresponding to an incremental volume of pumped fluid;
 - and pump control means, responsive to selection of a fluid delivery rate and responsive to the current sensing means, for operating the motor until the incremental rotation is sensed and for repeating the operation at variable time intervals selected in accordance with the selected fluid delivery rate.
- 45 8. A method for controlling the rate of fluid delivery in a medical fluid delivery system comprising a motor unit and a disposable fluid delivery set for mounting to the motor unit to form a pump, and wherein the motor unit is adapted to act on the fluid delivery set in repeated operating cycles to deliver a volume of fluid for each cycle and means are provided for detecting the completion of the operating cycle, which method comprises:
 - 50 operating the motor unit until the completion of the operating cycle is detected;
 - and repeating operation of the motor unit at time intervals selected in accordance with the desired rate of fluid delivery.
9. A method as claimed in claim 8 for further controlling the volume of fluid delivered comprising counting the number of the operating cycles and interrupting the repeating operation when a number of the
 - 55 operating cycles corresponding to a desired volume has been completed.
10. A method for controlling the rate of fluid delivery by a fluid delivery system comprising a motor unit and a disposable delivery set, wherein the pump unit includes a motor driven rotor for engaging a flexible pump tube on the delivery set to form a peristaltic pump and wherein at least one magnet is provided on

the rotor and a magnetic field detector is provided on the pump unit at a position adjacent the magnet in one rotational position of the rotor, which method comprises:
operating the pump unit to rotate the rotor until the magnetic field detector detects the magnet; and
repeating the operating step at a variable time interval selected according to the desired delivery rate of the
5 delivery system.

11. A method as claimed in claim 10 wherein there are provided a plurality of the magnets equally angularly spaced on the rotor and wherein the operating step comprises operating the pump unit until the magnetic field detector detects a selected number of the magnets.

12. A method as claimed in claim 10 or 11 wherein the operating step further includes the step of
10 inhibiting operation of the magnetic field detector for a selected time.

13. A method for controlling the rate of fluid delivery by a fluid delivery system comprising a motor unit and a disposable delivery set, wherein the pump unit includes a rotor driven by a DC motor for engaging a flexible pump tube on the delivery set to form a peristaltic pump, which method comprises:
detecting the AC component of the DC current and forming digital pulses representative thereof;
15 operating the pump unit to rotate the rotor until a selected number of the digital pulses occur; and
repeating the operating step at a time interval selected according to the desired delivery rate of the delivery system.

14. A method as claimed in claim 13 wherein the operating step further includes a step of inhibiting the AC detecting step for a selected time.

15. A method for detecting the absence of a pump tube on the rotor in a fluid delivery system comprising a motor unit and a disposable delivery set, wherein the pump unit includes a DC motor driven rotor for engaging a flexible pump tube on the delivery set to form a peristaltic pump, which method
20 comprises:

detecting the AC component of the current provided to the DC motor;
25 comparing the detected AC component of the current to a reference level thereby to detect current variation resulting from the presence of a pump tube on the rotor.

30

35

40

45

50

55

6-7593

APPENDIX I

Avocet 6805 Assembler v1.12, #01040 Chip=6305 12/18/87 16:08:57
K224/324 PUMP

```

1 $CHIP(6305)
2 $pagewidth=132
3 $PAGINATE
4 $TITLE(K224/324 PUMP)
5 $XREF
6 ;
7 ;
8 ; This file contains the software for the K224/K324 series of
9 ; Enteral Feeding Pumps. It is based on COMBO2.ASM a combined
10 ; F1500 and K224 program. See COMBO2.ASM for revision history.
11 ;
12 ;
13 ; FILE NAME= KPMP31.ASM
14 ;
15 ; DATE | REVISION | DESCRIPTION

```

```

16 ;-----
17 ;
18 ; 12/04/87 | 031.0 | DISPLAY LEADING ZERO'S IN DOSE MODE.
19 ; | | CHECK FOR DROPS CONTINUOUSLY IN RUN MODE
20 ; | | WHEN RATE IS < 100 ML/HR.
21 ;-----
22 ;
23 ; 11/24/87 | 030.0 | ADD NOPS TO IR TEST DURING TESTD.
24 ; | | CLEAR TESTD FLAG WHEN OFF MODE IS ENTERED TO
25 ; | | ALLOW PUMP TO DISABLE IF NECESSARY.
26 ; | | CLEAR VOL MODE FLAG WHEN LO BAT IS ENTERED
27 ; | | LOCK-OUT HOLD/START AND VOL BUTTONS WHEN
28 ; | | INC OR DEC BUTTONS ARE PRESSED.
29 ; | | FIX CLR VOL LED ERROR
30 ; | | CHECK DEC BUTTON BEFORE INC BUTTON.
31 ; | | ADD DELAY TO INC AND DEC BUTTON IN DOSE MODE TO
32 ; | | ALLOW TIME FOR DOSE DISPLAY.
33 ;-----
34 ;
35 ; 11/11/87 | 029.0 | CHANGE PUMP ID FROM PORTD(1), TO PORTB(6).
36 ;-----
37 ;

```



```

38 ; 10/23/87 | 028.0 | FIXED 24 HOUR TIME ERROR. PUMP WAS TIMING 21 HOURS
39 ; | | MODIFIED SO THAT PUMP DISABLES AFTER 24 HOURS
40 ; | | WHEN ON AC OR BATTERY POWER.
41 ; -----
42 ;
43 ; 10/16/87 | 027.0 | CHANGED SEICK ROUTINE TO TEST FOR SET WHENEVER
44 ; | | PUMP IS IN RUN MODE. OLD VERSION ONLY CHECKED
45 ; | | WHEN MOTOR IS RUNNING.
46 ; | | RATE CHECK FOR 324 IS DISABLED DURING TEST MODE.
47 ; | | DISABLE DOSE BUTTON WHILE INC OR DEC BUTTONS
48 ; | | ARE PRESSED.
49 ; -----
50 ;
51 ; 10/13/87 | 026.0 | ADD RECHECK OF PUMP ID DURING TESTD
52 ; | | ADDED DELAY FROM MOTOR TURN-ON TO DEAD BAT ROUTINE
53 ; | | CHANGED LOW BAT ROUTINE SO IT KILLS POWER IF
54 ; | | LO BAT OCCURS WHEN UNIT IS ALREADY OFF.
55 ; | | CHANGE POLARITY OF AC DETECT LINE. (REV 7 BOARD)
56 ; -----
57 ;
58 ; 10/13/87 | 025.0 | CORRECTED ERROR IN 1/2 SEC TIMER..
59 ; | | ADDED DOSE DEL CHECK PRIOR TO RUN MODE

```

5-7-83
1333

```

60 ;-----
61 ;
62 ; 10/12/87 | 024.0 | ADD NOPS IN DROPCK BETWEEN IR TURN-ON AND
63 ; | | TEST TO ALLOW FOR SETTLING.
64 ; | | SHORTED IR ON TIME.
65 ; | | MOVE DROPCK FROM 8 MS TO 2 MS TO IMPROVE
66 ; | | DROP DETECTION.
67 ;-----
68 ;
69 ; 10/10/87 | 023.0 | DEL STOP INSTR. FROM EACH DISABLE ROUTINE.
70 ; | | BLANK ALL DIGITS IN EACH DISABLE ROUTINE.
71 ;-----
72 ;
73 ; 10/07/87 | 022.0 | THE RESET ROUTINE HAS BEEN MODIFIED TO KEEP
74 ; | | THE PUMP OFF. TURN-ON IS NOW ACCOMPLISHED IN
75 ; | | THE INT ROUTINE. CHANGE WAS MADE TO PREVENT STRAY
76 ; | | RESEST SIGNALS FROM WATCHDOG FROM ACCIDENTLY
77 ; | | TURNING PUMP ON DURING DISABLE SEQUENCE.
78 ;-----
79 ; 9/29/87 | 021.0 | REVISED THE DISABLE ROUTINES TO DISABLE INT
80 ; | | AND LOOP UNTIL POWER DISSIPATES.
81 ;-----

```

```

82 ;
83 ; 9/22/87 | 020.0 | TEST ROUTINE HAS BEEN ADDED. IT IS INITIATED BY
84 ; | | PRESSING THE OFF AND HOLD BUTTONS DOWN AT THE
85 ; | | SAME TIME FOR APPROX. 3 SEC WHEN PUMP IS OFF.
86 ;
87 ; 9/22/87 | 019.0 | MOTOR CONTROL SIGNAL POLARITY HAS BEEN REVERSED
88 ; | | AND TEST MODE WAS CHANGED SO THAT CLR VOL LED
89 ; | | IS TURNED OFF WHEN BUZZER IS.
90 ;
91 ; 9/08/87 | 018.0 | THE INCREMENT AND DECREMENT KEYS HAVE BEEN
92 ; | | REVERSED TO COMPENSATE FOR PC BOARD R4.
93 ; -----
94 ;
95 ; 9/03/87 | 017.0 | THIS PROGRAM HAS BEEN MODIFIED TO RUN ON A 2 MHZ
96 ; | | CLOCK.
97 ; | | THE LOW BAT TIME HAS BEEN SET TO 15 MIN.
98 ; -----
99 ; -----
100 ;
101 ; 8/28/87 | 016.0 | THE CLEAR DOSE FEATURE HAS BEEN DELETED
102 ; -----
103 ;

```

```

104 ; 8/28/87 | 015.0 | THE HOLD ROUTINE HAS BEEN MODIFIED TO OCCUR
105 ; | | 2 1/2 MIN FROM LAST BUTTON PRESSED.
106 ; | | A WAIT STATEMENT HAS BEEN ADDED TO THE MAIN LOOP
107 ; | | AND OFF LOOP ROUTINES.
108 ; -----
109 ; 8/21/87 | 014.0 | THIS VERSION HAS BEEN MODIFIED TO ACCOMADATE
110 ; | | THE NEW ENCODED ROTARY SWITCH.
111 ; | | THE HOLD AND INCREMENT BUTTIONS HAVE BEEN SWITCHED
112 ; | | TO SIMPLIFY THE RATE KNOB DECODE OPERATION.
113 ; -----
114 ;
115 ; 8/12/87 | 013.0 | THIS VERSION HAS BEEN REDUCED BY CONVERTING
116 ; | | THE RATE INCREMENT ROUTINES TO BINARY WITH
117 ; | | DECIMAL ADJUST ADDING.
118 ; -----
119 ;
120 ; 8/05/87 | 012.0 | BLINK SUBROUTINE WAS MODIFIED SO THAT ONLY 3
121 ; | | DIGITS ARE ACTIVE WHEN PUMP IS 224.
122 ; -----
123 ;
124 ; 7/30/87 | 011.0 | MODIFIED SO THAT VOLUME AND DOSE CAN BE CLEARED
125 ; | | WITH CLR V/D BUTTON WHILE PUMP IS RUNNING.


```

```

126 ;      |      |      | THE DISPLAY HAS BEEN CHANGED TO EYE ERR IF
127 ;      |      |      | THE IR SENSOR IS BLOCKED DURING TEST SEQUENCE.
128 ;-----
129 ;
130 ; 7/25/87 | 010.0 | MODIFY DOSE DEL ERROR, SO THAT PRESSING HOLD
131 ;      |      |      | BUTTON ONCE DISABLED ALARM AND MESSAGE. THIS
132 ;      |      |      | WILL ALLOW THE USER TO UPDATE DOSE OR VOLUME
133 ;      |      |      | WITH THE NEXT KEY STROKE.
134 ;-----
135 ;
136 ; 7/21/87 | 009.0 | modified to display DP scroll when motor on
137 ;      |      |      | in LOW BAT state.
138 ;-----
139 ;
140 ; 7/16/87 | 008.0 | MODIFY HOLD BUTTON OPERATION FOR LO BAT STATE.
141 ;      |      |      | WHEN HOLD PRESSED, ALARM STOPS AND MOTOR STOPS.
142 ;-----
143 ;
144 ; 7/13/87 | 007.0 | ADD NO SET FEATURE.
145 ;      |      |      | MODIFY VOLUME DISPLAY. ADD SHORTED MOTOR SENSOR
146 ;      |      |      | DETECTION SOFTWARE. MODIFY 24 HR TIME ROUTINE
147 ;      |      |      | TO CLEAR VOL DELIVERED, BUT NOT TO KILL POWER.

```

```
148 ;-----
149 ; 7/10/87 | 006.0 | FIX HIGH RATE ERROR BUG. MODIFY LOW BAT TO
150 ; | | CONTINUE TO RUN MOTOR WITH ALARM.
151 ;-----
152 ;
153 ; 7/08/87 | 005.0 | TRY NEW INCREMENT METHOD FOR DOSE.
154 ; | | SAVES APPROX. 40 BYTES. ADDED LOW BAT TIMER,
155 ; | | AND DEAD BAT SIGNAL.
156 ;-----
157 ;
158 ; 7/07/87 | 004.0 | 324 RATE INCREMENT CHANGED TO 5 ML FOR RATES
159 ; | | GREATER THEN 50.
160 ;-----
161 ; | |
162 ; 6/29/87 | 003.0 | DOSE FEATURE IS ADDED. CURRENT SENSE SOFTWARE FOR
163 ; | | SET DETECTION DELETED.
164 ; | |
165 ;-----
166 ; 6/25/87 | 002.0 | RATE INITIALIZED TO 0 ON POWER UP AND
167 ; | | TURN-ON. WRAP AROUND OF DISPLAY REMOVED.
168 ;
169 ;-----
```

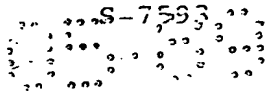

| | | | | |
|-----|---|---------|-------|------------------------------|
| 170 | ; | ; | ; | ; |
| 171 | ; | 6/01/87 | 001.0 | FIRST PUMP VERSION SEPARATED |
| 172 | ; | ; | ; | FROM COMBO2.ASM |
| 173 | ; | ; | ; | |
| 174 | | \$EJECT | | |

```

175 ;
176 ;
177 ;
178 ; PORT ASSIGNMENTS A ALL OUTPUTS
179 ;
180 ;
181 ; -----
182 ; 7 6 5 4 3 2 1 0
183 ; DP /SEG g /SEG f /SEG e /SEG d /SEG c /SEG b /SEG a
184 ; -----
185 ;
186 ;
187 ; B ALL OUTPUTS
188 ;
189 ;
190 ; -----
191 ; 7 6 5 4 3 2 1 0
192 ; BUZZER VOL CLR 5VA POWER /DIGIT 4 /DIGIT 3 /DIGIT 2 /DIGIT 1
193 ; EMPID 1=224 1=ON ON/OFF (324)
194 ; -----
195 ;
196 ;

```



```

197 ;
198 ;
199 ;
200 ;
201 ; 7 6 5 4 3 2 1 0
202 ; F1500/ AC/DC /DROP /IO BAT /SET /OFF /B /A /A
203 ; PUMP
204 ;
205 ;
206 ;
207 ;
208 ;
209 ;
210 ;
211 ; 7 6 5 4 3 2 1 0
212 ; NC IR-LED NC CLEAR VOLUME /HOLD /DOSE 1/3
213 ; ON/OFF REV
214 ;
215 ;
216 ;
217 ;
218 $EJECT

```

C ALL INPUTS

D 0,1,2,3,4, INPUTS
5,6 OUTPUTS

```

219      DEFSEG ABSLSEG,ABSOLUTE
220      SEG ABSLSEG
221      ;
222      ;
223      ;***** SET UP RAM/ROM AND PORT ADDRESSES *****
224      ;
225      USRRAM EQU $40      ;STARTING ADDRESS OF USER RAM.
226      USRRROM EQU $1000   ;STARTING ADDRESS OF USER ROM.
227      EMPTST EQU $1C00    ;STARTING ADDRESS OF PUMP TEST SOFTWARE.
228      ;
229      PORTA EQU $00        ;PORT A DATA.
230      PORTB EQU $01        ;PORT B DATA.
231      PORTC EQU $02        ;PORT C DATA.
232      PORTD EQU $03        ;PORT D DATA.
233      ;
234      DDRA EQU $04         ;PORT A DATA DIRECTION REGISTER.
235      DDRB EQU $05         ;PORT B " " " "
236      DDRC EQU $06         ;PORT C " " " "
237      DDRD EQU $07         ;PORT D " " " "
238      TDR EQU $08          ;TIMER DATA REGISTER.
239      TCR EQU $09          ;TIMER CONTROL REGISTER.
240      MISC EQU $0A         ;MISCELLANEOUS REGISTER

```

```

241 ;
242 ; *****
243 ;          EQUATES
244 ; *****
245 ;
246      TABLE1 EQU $101C
247      TABLE5 EQU $1000
248 ;
249      TIMEST EQU $02
250 ;
251      T8MS EQU $03
252      T33MS EQU $03
253      T1S EQU $0F          ;1/2 SEC TIME CONSTANT
254      T6M EQU $02          ;(02 + 1) = 03 x 2.0 = 6 MINUTES
255      T142S EQU $EF         ;(EF + 1) = 240 x 0.5 = 2 MINUTES
256 ;
257 ;
258 ;***** INTERRUPT/RESET VECTOR TABLE *****
259 ;
260      ORG $1FF4
261 ;
262      1FF4 102D            DW RESET ;SERIAL INTERRUPT TIMER 2

```

| | | | | |
|-----------|-----|----|--------|--------------------------------------|
| 1FF6 1089 | 263 | DW | TIMINP | ;TIMER INTERRUPT VECTOR (WAIT STATE) |
| 1FF8 1089 | 264 | DW | TIMINP | ;TIMER INTERRUPT VECTOR |
| 1FFA 103E | 265 | DW | ONPOW | ;EXTERNAL INTERRUPT VECTOR |
| 1FFC 103E | 266 | DW | ONPOW | ;SOFTWARE INTERRUPT VECTOR |
| 1FFE 102D | 267 | DW | RESET | ;RESET VECTOR |
| | 268 | | | ; |
| | 269 | | | ; |
| | 270 | | | \$EJECT |

```

271 ;
272 ;***** DEFINE VARIABLES IN USER RAM AREA *****
273 ;
274 ;
275 =0040          ORG      USRRAM
276 ;
277 TS            DB      0      ;READING OFF BRIDGE BEFORE RISE BEGINS (AMBIENT)
278 T5            DB      0      ;READING FROM AMP 5 SECOND AFTER BRIDGE RISE.
279 T10           DB      0      ;      "      10
280 T15           DB      0      ;      "      15
281 T20           DB      0      ;      "      20
282 T25           DB      0      ;      "      25
283 ;
284 ALGO          DB      0
285 ;
286 TIMLO         DB      0      ;LOWER BYTE OF 5 MINUTE DOWN COUNTER.
287 TIMHI         DB      0      ;UPPER " " " "
288 ;
289 BATIM1         DB      0      ;15 MIN BATTERY TIMERS
290 BATIM2         DB      0
291 ;
292 HLDTIM1        DB      0      ;HOLD TIME COUNTERS: 2 1/2 MIN

```

```

004C 00      293 HLDIM2 DB 0
                294 ;
004D 00      295 SAMPL1 DB 0 ;STORAGE FOR THE LAST 5 A-TO-D READINGS.
004E 00      296 SAMPL2 DB 0 ;(SAMPLE5 IS THE MOST RECENT).
004F 00      297 SAMPL3 DB 0
0050 00      298 SAMPL4 DB 0
                299 ;
0051 00      300 DSPDG1 DB 0 ;7-SEGMENT DATA FOR LSD. (BIT 7 IS BACKPLANE CLOCK)
0052 00      301 DSPDG2 DB 0 ;" " " " DIG 2.
0053 00      302 DSPDG3 DB 0 ;" " " " DIG 3. (BIT 7 IS POWER CONTROL)
0054 00      303 DSPDG4 DB 0 ;" " " " MSD. (BITS 7-3 ARE A/D CONTROL)
                304 ;
                305 ;
                306 ;
0055 00      307 DS1 DB 0 ;CONTAINS FIRST TWO DIGITS OF DOSE
0056 00      308 DS2 DB 0 ;CONTAINS LAST TWO DIGITS OF DOSE
                309 ;
                310 ; DS1 AND DS2 ARE TRANSFERRED TO SAMPLE1-SAMPLE4 BY BCDEXP
                311 ; ROUTINE.
                312 ;
0057 00      313 DOSECT DB 0 ;COUNTER FOR TIME IN DOSE MODE
                314 ;

```

```

315 ; VS1 AND VS2 ARE TRANSFERRED TO SAMPL1-SAMPL4 BY BCDEXP
316 ;
317 DECML1 DB 0 ;FRACTION OF 1ML IN VOLUME
318 VS1 DB 0 ;CONTAINS FIRST TO DIGITS OF VOLUME
319 VS2 DB 0 ;CONTAINS LAST TO DIGITS OF VOLUME
320 ;
321 REG1 DB 0 ;TEMPORARY STORAGE
322 ;
323 REG2 DB 0 ; "
324 ;
325 DRCNT1 DB 0 ;DROP COUNTERS
326 DRCNT2 DB 0
327
328 REG3 DB 0 ; "
329 REG4 DB 0 ; "
330 REG5 DB 0 ; "
331 ;
332 ;
333 ;
334 ;
0058 00
0059 00
005A 00
005B 00
005C 00
005D 00
005E 00
005F 00
0060 00
0061 00

```

```

0062 00      COUNT1 DB 0      ;GENERAL PURPOSE COUNTER REGISTER.
0063 00      COUNT2 DB 0      ;      "      "      "
0064 00      COUNT3 DB 0      ;      "      "      "
338 ;
339 ;***** REGISTERS AND VARIABLES USED DURING PREDICTION *****
340 ;
341 ;
342 XH DB 0      ;ARGUMENT STORAGE FOR 16-BIT MATH ROUTINES.
343 XL DB 0      ;      "
344 QH DB 0      ;      "
345 QL DB 0      ;      "
346 PH DB 0      ;      "
347 PL DB 0      ;      "
348 ;
349 TEMPA DB 0
350 TEMPX DB 0
351 ;
352 PSTAT1 DB 0      ;FIRST PUMP STATUS BYTE
353 ;      0=HOLD      1=Unit on hold
354 ;      1=MOTOR      1=Motor on
355 ;      2=ON      1=Unit on
356 ;      3=START      1=Unit in start mode

```



```

357 ;
358 ;
359 ;
360 ;
361 ;
362 ;
363 PSTAT2 DB 0 ;SECOND PUMP STATUS BYTE
364 ;
365 LOWBAT EQU PSTAT2 0=BATTERY 1=Low battery
366 ACON EQU PSTAT2 1=AC 1=AC ON
367 NDROP EQU PSTAT2 2=FLOW 1=Drop error
368 HOLDER EQU PSTAT2 3=HOLD 1=Hold error
369 NOSET EQU PSTAT2 4=SET 1=No set error
370 ALRM EQU PSTAT2 5=ALARM 1=Audio alarm on
371 VOIM EQU PSTAT2 6=VOLUME 1=Volume mode
372 ERRON EQU PSTAT2 7=ERROR 1=Error mode
373 ;
374 TSTERR DB 0 ;B0 IS TEST ERROR FLAG
375 ; SET IF EYES BLOCKED DURING TEST
376 ;
377 DOSER DB 0 ;DOSE DELIVTRED FLAG B0=1: DOSE DEL
378 DFLAG DB 0 ;DOSE MODE FLAG


```

```

379 ; B0= DOSE 1=DOSE MODE ON
380 ; B1= DTEST 1=DOSE BUTTON PRESSED
381 ; LAST CYCLE, USED FOR
382 ; DEBOUNCE.
383 ; B2= ZERO 1=DISPLAY ZERO
384 ;
0072 00 LTEST DB 0 ;MAGNET LOW COUNTER, USED TO CHECK IF SENSOR
385 ; IS STUCK LOW. INCREMENTED AS LONG REV SIGNAL
386 ; IS LOW. CLEARED WHEN SIGNAL GOES HIGH.
0073 00 HTEST DB 0 ;B0 - TEST FLAG FOR HOLD DEBOUNCE
387 ;
388 ; B1 - TEST FLAG FOR CIR VOL DEBOUNCE
389 ;
0074 00 ZEROST DB 0 ;ZERO STATUS FLAG
390 ;
391 ; B5 = 1 : DISPLAY IS ZERO
392 ; B6 = 1 : MSD IS ZERO
393 ; B7 = 1 : MSD-1 IS ZERO
394 ;
0075 00 VFLAG DB 0 ;VOLUME DISPLAY TEST FLAG. SET WHEN 'VOL'
395 ; IS DISPLAYED.
396 ;
0076 00 CWF DB 0 ;CLOCKWISE FLAG
0077 00 CCWF DB 0 ;COUNTER CLOCKWISE FLAG
0078 00 SPEED1 DB 0 ;ROLL UP/DOWN SPEED
0079 00 VOLTIM DB 0 ;VOLUME MODE ON TIMER

```

3NSDOCID: <EP_____0327209A2_1_>

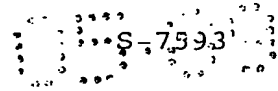


```

419 ; *****
420 ;          CONSTANTS IN ROM
421 ; *****
422 ;
      LO DB $F8,$9C,$FF ;DISPLAY LO ON LEDS
      BATT DB $98,$81,$B8 ;DISPLAY BAT ON LEDS
      FLO DB $B1,$F8,$C0 ;DISPLAY FLO ON LEDS
      ERR DB $B0,$BD,$BD ;DISPLAY ERR ON LEDS
      NO DB $9D,$9C,$FF ;DISPLAY NO ON LEDS
      SET DB $92,$B0,$B8 ;DISPLAY SET ON LEDS
      HLD DB $89,$F8,$8C ;DISPLAY HLD ON LEDS

1000 F8 9C FF
1003 98 81 B8
1006 B1 F8 C0
1009 B0 BD BD
100C 9D 9C FF
100F 92 B0 B8
1012 89 F8 8C

1015 8C 9C 92 B0 DOSE DB $8C,$9C,$92,$B0 ;DISPLAY DOSE ON LEDS
1019 8C B0 F8 DEL DB $8C,$B0,$F8 ;DISPLAY OUT ON LEDS
101C 92 8A 92 SYS DB $92,$8A,$92 ;DISPLAY SYS ON LEDS
101F C0 CF A4 86 DIGITS DB $C0,$CF,$A4,$86,
1023 8B $8B ;DISPLAY NUMBERS 0-9 ON LED
1024 92 90 C7 80 DB $92,$90,$C7,$80,
1028 83 $83
435 ;
436 DNUM DB $F7,$FB,$FD,$FE ;DIGIT OUTPUT

```

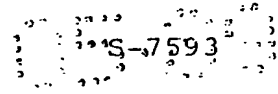

437 ;
438 ;*****
439 ;
440 \$EJECT

3-7593

```

441 ;
442 ;-----
443 ;          START FROM RESET
444 ;-----
445 ;
446 ;          (THIS IS THE POWER-ON RESET ENTRY POINT)
447 ;
448 ;*****
449 ;
450 ;
451 102D AE 40      RESET LDX    #$40      ;FIRST, CLEAR RAM.
452 102F 4F          CLR A
453 1030 F7      CLEAR STA    ,X
454 1031 5C          INCX
455 1032 A3 FF      CPX    #$FF
456 1034 26 FA      BNE    CLEAR
457
458 ;
459 ;          Set up the ports
460 ;
461 1036 CD 1106      JSR    TIMEON      ;GO START TIMER COUNTER
462 1039 1A 0A      BSET 5,MISC      ;SET INT SENSE FOR EDGE AND LEVEL

```

```

463      ;
464      ;      Wait for timer interrupt
465      ;
466      LOOPW  WAIT
103B 8F
103C 20 FD      BRA  LOOPW
467
468      ;
469      $EJECT

```

```

470 ;
471 ;
472 ;
473 ;
474 ;
475 ;
476 ;
477 ;
478 ONPOW BC1R 5,MISC ;RETURN INT SENSE TO EDGE LEVEL
479 BRSET 7,ERRON,NOON ;IF IN ERROR MODE DONT DO
480 BRSET 2,PSTAT1,NOON ;DONT DO IF ON
481 ;
482 LDA #19 ;INITIALIZE PUMP STAT REGISTER
483 STA PSTAT1
484 ;
485 ;
486 ;
487 BSET 7,PSTAT1 ;PRESET ID FLAG FOR 324
488 LDA #FF ;PRELOAD ACCA TO SET PORTB ALL OUTPUT
489 BRCLR 6,PORTB,CONON ;IF PORTB(6)=0, PMP IS 324 SO CONTINUE
490 BC1R 7,PSTAT1 ;ELSE, PMP IS 224 SO CLEAR ID FLAG
491 LDA #BF ; AND SET ACCA TO MAKE PORTB(6) INPUT

```

EXTERNAL INTERRUPT HANDLER

```

103E 1B 0A
1040 0E 6E 45
1043 04 6D 42
1046 A6 19
1048 B7 6D
104A 1E 6D
104C A6 FF
104E 0D 01 04
1051 1F 6D
1053 A6 BF

```



```

492      ;
493      CONON STA      DDRB      ;CONFIGURE PORTB TO BF OR FF
494      LDA      #$FF
495      STA      DDRA      ;CONFIGURE PORTA TO ALL OUTPUTS
496      STA      PORTA      ;MAKE PORTA ALL 1'S
497      LDA      #$1F      ;
498      STA      PORTB      ; PORTB=0001 1111
499      LDA      #$60      ;MAKE PORT D 0-4 INPUTS 5,6 OUTPUTS
500      STA      DDRD
501      CLRA
502      STA      PORTD
503      STA      DDRC
504      ;
505
506      ;      Initialize display pointers
507      ;
508      ;SET UP INITIAL DISPLAY
509      ;POINTERS SO DISPLAY
510      ;READS 0
511      ;WHEN THE DISPLAY IS
512      ;ACTIVATED
513      ;

```

ISDOCID: <EP_____0327209A2 I_>

| | | | |
|---------|-----|---------|------|
| 1088 80 | 536 | NOON | RITI |
| | 537 | ; | |
| | 538 | ; | |
| | 539 | \$EJECT | |

VSDOCID: <EP_____0327209A2_I_>

```
1095 CD 12DA      562      JSR  INRATE1      ;GO UPDATE RATE FOR 224
1098 CD 1185      563      JSR  DROPCK      ;GO CHECK FOR DROP
                    564      ;
109B B6 40        565      LDA  TS          ;TEST IF TIME TO DO NEXT TESTS
109D 3C 40        566      INC  TS
109F A1 03        567      CMP  #T8MS
10A1 27 01        568      BEQ  TIME2
10A3 80           569      RTI
                    570      ;
                    571      ; ***** 8 MS LOOP *****
                    572      ;
10A4 3F 40        573      TIME2 CLR TS
                    574      ;
10A6 CD 11CF      575      JSR  RPMCK      ;GO CHECK MOTOR FOR ONE RPM
10A9 CD 1290      576      JSR  MOTCK      ;GO CHECK MOTOR TIME TO TURN ON
10AC CD 12BD      577      JSR  SETCK      ;GO CHECK FOR SET CURRENT
10AF CD 140E      578      JSR  OFFCK      ;GO SEE IF OFF PUSHED
                    579      ;
10B2 B6 41        580      LDA  T5          ;GET STATUS OF SECOND TIMER
10B4 3C 41        581      INC  T5          ;INCREMENT FOR NEXT TIME
10B6 A1 03        582      CMP  #T33MS      ;DOES IT EQUAL 32 MS
10B8 27 01        583      BEQ  TIME3      ;YES SO GO TEST
```

6-17593

```

10BA 80      584      RTI      ;RETURN TO PROGRAM
              585      ;
              586      ; ***** 33 MSEC CHECK LOOP *****
              587      ;
10BB 3F 41    588      TIME3 CLR T5      ;CLEAR T5 FOR NEXT TIME THRU
              589      ;
10BD CD 1458  590      JSR TESTD      ;DO TEST DISPLAY
10C0 CD 14EC  591      JSR HOLDCK      ;GO TEST IF HOLD PRESSED
10C3 CD 15A9  592      JSR DOSECK      ;GO CHECK DOSE BUTTON
10C6 CD 15EE  593      JSR VOLDIS      ;DISPLAY VOLUME CHECK MODEL 324
10C9 CD 1640  594      JSR CKCLR      ;GO CHECK VOLUME CLEAR
10CC CD 1687  595      JSR INRATE2      ;INPUT ROUTINE FOR MODEL 324
              596      ;
10CF B6 42    597      LDA T10      ;GET TIMER3 STATUS
10D1 3C 42    598      INC T10      ;INCREMENT FOR NEXT TIME
10D3 A1 0F    599      CMP #T1S      ;DOES IT EQUAL 1S
10D5 27 01    600      BEQ TIME4      ;YES SO GO TIME4
10D7 80      601      RTI      ;RETURN TO PROGRAM
              602      ;
              603      ; ***** 1/2 SEC BAT CHECK *****
              604      ;
10D8 3F 42    605      TIME4 CLR T10      ;RESET TIMER3 FOR NEXT TIME

```

```

606      ;
        10DA CD 17DD      JSR  ADCD      ;GO TEST AC OR DC
607      ;
        10DD CD 17E6      JSR  BATCK     ;GO DO BATTERY CHECK
608      ;
        10E0 CD 1836      JSR  DBATCK    ;GO CHECK FOR DEAD BATTERY
609      ;
        10E3 CD 184E      JSR  ALARM     ;GO TEST ALARM
610      ;
        10E6 CD 1860      JSR  BLINK     ;GO BLINK DISPLAY
611      ;
        10E9 CD 1909      JSR  HILDER    ;go test if on hold 142sec then error
612      ;
        613      ;
        614      NO5      LDA  T15      ;GET TIMER4 STATUS
615      ;
        10EC B6 43      INC  T15      ;ADD ONE FOR NEXT TIME
616      ;
        10EE 3C 43      CMP  #T142S    ;DOES IT EQUAL 142 SEC
617      ;
        10F0 A1 EF      BEQ  TIMES5     ;YES SO GO TIME 5
618      ;
        10F2 27 01      RTI            ;RETURN TO PROGRAM
619      ;
        620      ;
        621      ;
        622      TIMES5 CLR  T15      ;CLEAR FOR NEXT TIME
623      ;
        624      ;
        625      LDA  T20
626      ;
        10F7 B6 44      INC  T20
627      ;
        10F9 3C 44      CMP  #T6M
        10FB A1 02

```

***** 142 SEC LOOP *****

↓SDOCID: <EP_____0327209A2_I_>


37593

```

641 ;
642 ; SUBROUTINES
643 ;
644 ;
645 ;
646 ;
647 ; INITIALIZE COUNTER TIMER
648 ;-----
649 ;
650 ;
651 ;
652 TIMEON LDA #TIMEST ;SET COUNTER TO 4
653 STA TCR ;ENABLE TIMER INT
654 CLI ;ENABLE INTERRUPTS
655 RTS ;RETURN;
656 ;
657 ;
658 ;
659 ; WATCH DOG TIMER
660 ;-----
661 ; Output clock to hold off reset
662 ;

```

1106 A6 02
1108 B7 09
110A 9A
110B 81

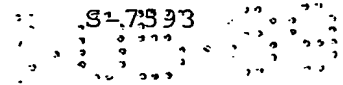

| | | | | |
|---------------|-----|--------|-----------------------|-------------------------|
| 110C 0B 03 03 | 663 | WATCHD | BRCIR 5, PORTD, SETTR | ;IF CLOCK LOW MAKE HIGH |
| 110F 1B 03 | 664 | | BCLR 5, PORTD | |
| 1111 81 | 665 | | RTS | |
| | 666 | | ; | |
| 1112 1A 03 | 667 | | SETTR BSET 5, PORTD | |
| 1114 81 | 668 | | RTS | |
| | 669 | | ; | |
| | 670 | | \$EJECT | |

```

671 ;
672 ;
673 ;
674 ;
675 ;
676 ;
677 ;
678 ;
679 ;
680 ;
681 ;
682 ;
683 ;
684 ;
685 ;
686 ;
687 ;
688 ;
689 ;
690 ;
691 ;
692 ;

                                LED MUX ROUTINE

                                This routine is called every 4 ms
                                It handles 0 suppression and outputs 7 seg data
                                VARIABLES USED: COUNT3=This is the digit number and is changed every
                                                time this routine is called.
                                                SAMPLE1=location of Current value of digit 1.
                                                add count3 and use as pointer.
                                                7,pstat2 (error)=test if in error mode for blinking.
                                                4,pstat1 (display)=test if display blanked (blinking
                                                mode).
                                                ZEROSET=bit 6 set=zero in msd
                                                bit 7 set=zero in next digit

                                -----

                                LED    BRCLR    2,PSTAT1,NOD4    ;IF UNIT OFF DONT DO MUX
                                LDA    #$0F                    ;Turn all digits off
                                ORA    PORTB

```

1115 05 6D 36
1118 A6 0F
111A BA 01

```

111C B7 01      693      STA      PORTB
111E BE 64      694      LDX      COUNT3      ;GET CURRENT DIGIT NO.
1120 A3 03      695      CPX      #$03      ;HAVE WE DONE 4
1122 26 04      696      BNE      GO      ;IF NOT GO ON
1124 3F 64      697      CLR      COUNT3      ;IF SO RESET ALL VARIABLES
1126 20 02      698      BRA      GO2
699      ;
1128 3C 64      700      GO      INC      COUNT3
701      ;
112A 0C 6E 0E   702      BRSET 6,VOLM,NOTZRO      ;IF IN VOL NO ZERO TEST
112D 00 71 0B   703      BRSET 0,DFLAG,NOTZRO      ;IF IN DOSE NO ZERO TEST
1130 A3 03      704      CPX      #$03      ;IF DIGIT 1 DONT TEST FOR 0
1132 27 07      705      BEQ      NOTZRO
1134 CD 1175     706      JSR      TEST0      ;GO TEST
1137 24 02      707      BCC      NOTZRO      ;IF NO CARRY IT CANT BE 0
1139 20 14      708      BRA      ZERO
709      ;
710      ;
113B E6 4D      711      NOTZRO      LDA      SAMPLE1,X
113D 0E 6E 05   712      BRSET 7,ERRON,OUT1      ;IF error mode dont blank
1140 08 6D 02   713      BRSET 4,PSTAT1,OUT1      ;if display on flag set dont blank
1143 4F         714      CIRA      ;blank display if display flag is low

```

```

1144 43      715      COMA      ;SET PORTA=FF ALL OFF
              716      ;
1145 B7 00   717      OUT1 STA  PORTA
1147 D6 1029 718      LDA  DNUM,X      ;GET DIGIT TO OUTPUT
114A B4 01   719      AND  PORTB      ;ENABLE
114C B7 01   720      STA  PORTB
114E 81      721      NOD4 RTS
              722      ;
              723      ;
114F A3 00   724      ZERO CPX  #$00      ;IS THIS FIRST TIME THROUGH (MSD)
1151 27 17   725      BEQ  ZERO1      ;IF SO BLANK AND SET FLAG IF ZERO
1153 A3 01   726      CPX  #$01      ;IS THIS SECOND TIME THRU
1155 27 17   727      BEQ  ZERO2      ;IF SO CHECK FIRST FLAG
1157 0D 74 E1 728      BRCLR 6,ZEROST,NOTZRO ;ON THIRD TIME CHECK FIRST TWO
115A 0F 74 DE 729      BRCLR 7,ZEROST,NOTZRO ;IF EITHER NOT SET DONT BLANK
              730      ;
115D E6 4D   731      ZEROT LDA  SAMPL1,X  ;GET DIGIT DATA
115F 49      732      ROLA      ;CHECK IF DP ON
1160 24 04   733      BCC  ZERODP      ;IF SO LEAVE ON
1162 A6 FF   734      LDA  #$FF      ;IF NOT BLANK
1164 20 DF   735      BRA  OUT1      ;GO OUTPUT IT
1166 A6 7F   736      ZERODP LDA  #$7F      ;LEAVE DP ON

```

```

1168 20 DB          BRA OUT1
737
738 ;
739 ;
116A 1C 74          ZERO1 BSET 6,ZEROST ;SET ZERO FLAG
116C 20 EF          BRA ZEROT ;OUTPUT BLANK
740
116E 0D 74 CA          ZERO2 BRCIR 6,ZEROST,NOTZRO ;CHECK FIRST DIGIT FLAG
741
1171 1E 74          BSET 7,ZEROST ;IF SET THEN BLANK 2ND DIGIT
742
1173 20 E8          BRA ZEROT
743
744
745 ;
1175 E6 4D          TEST0 LDA SAMPL1,X ;test if digit is 0
746
1177 98          CLC
747
1178 49          ROLA
748
1179 A1 80          CMP #80
749
117B 27 04          BEQ SETCAR
750
117D 98          CLC ;clr carry if not zero
751
117E E6 4D          LDA SAMPL1,X
752
1180 81          RIS
753
1181 99          SETCAR SEC ;set carry if it is
754
1182 E6 4D          LDA SAMPL1,X
755
1184 81          RIS
756
757 ;
758 $EJECT

```

```

759 ;
760 ; DROP CHECK
761 ;
762 ; This routine is called every 8ms
763 ; After a delay of .4 sec to clear all transients,
764 ; it turns on the IR source and looks at the photo transistor.
765 ; If the line is low a drop is there the counter is cleared.
766 ; If there is no drop the counter is incremented and if it reaches
767 ; a count of 95 or 95 X .008=.76 sec then the second counter is incremented
768 ; This counter is tested in the RPM routine if it gets to 2 or 1.76s then
769 ; it will alarm.
770 ;
771 ; -----
772 ; Conditions
773 ;
1185 05 6D 43 DROPCK BRCLR 2,PSTAT1,NODRP ;DONT DO IF UNIT OFF
1188 00 6D 40 BRSET 0,PSTAT1,NODRP ;DONT DO IF PUMP NOT IN RUN MODE
776 ;
777 ; CHECK RATE
778 ;
118B B6 53 LDA DSPDG3 ;IF RATE < 100, CHECK FOR DROPS
118D A1 01 CMP #501

```

```

118F 25 03      781      BLO DELCK
1191 03 6D 37      782      BRCLR 1,PSTAT1,NODRP      ELSE, CHECK DROPS WHEN MOTOR ON ONLY
                        783      ;
                        784      ;          DELAY CHECK
                        785      ;
1194 B6 5B      786      DELCK LDA REG1          ;CHECK REG1
1196 A1 7C      787      CMP #124          ;IF 0.248 SEC ELAPSED SINCE MOTOR ON
1198 27 04      788      BEQ DROP1          ;DO DROP TEST AND SKIP INCREMENT.
119A 4C      789      INCA          ;ELSE,
119B B7 5B      790      STA REG1          ;INCREMENT AND RETURN.
119D 81      791      RTS
                        792      ;
                        793      ;          Test for drop
                        794      ;
119E 1C 03      795      DROP1 BSET 6,PORTD          ;TURN IR SOURCE ON
11A0 9D      796      NOP          ;ADD SOME DELAY TO ALLOW SETTLING
11A1 9D      797      NOP
11A2 9D      798      NOP
11A3 9D      799      NOP
                        800      ;
11A4 0A 02 1C      801      BRSET 5,PORTC,CLERK          ;CHECK FOR LOW (DROP)
11A7 1D 03      802      BCLR 6,PORTD          ;TURN OFF IR SENSOR

```


5-7593

```

11A9 10 81      803      BSET 0,COUNT4      ;WE HAVE DROP SO SET FLAG
                        804      ;
11AB 3C 5D      805      INC DRCNT1      ;INCREMENT DROP COUNTERS
11AD 26 02      806      BNE CONDRP
11AF 3C 5E      807      INC DRCNT2
11B1 B6 5E      808      CONDRP LDA DRCNT2      ;CHECK IF DROP CNT = 0.76 SEC
11B3 A1 01      809      CMP #01
11B5 25 14      810      BLO NODRP      ;IF DRCNT2 < 1 END TEST
11B7 B6 5D      811      LDA DRCNT1      ;ELSE CHECK DROP CNT 1
11B9 A1 7C      812      CMP #7C
11BB 27 02      813      BEQ NEXTBY      ;IF DRCNT1 >= 7C, INCREMENT REG3
11BD 20 0C      814      BRA NODRP      ;ELSE, EXIT DROPCK
                        815      ;
11BF 3C 5F      816      NEXTBY INC REG3
11C1 20 08      817      BRA NODRP      ;REG3 IS TESTED IN THE RPM ROUTINE
                        818      ; IF 2 IT ERRORS
                        819      ;
                        820      ; No drop
                        821      ;
11C3 1D 03      822      CLERK BCLR 6,PORTD
11C5 3F 5D      823      CLR DRCNT1      ;CLEAR THE COUNTERS FOR BLOCKED EYES
11C7 3F 5E      824      CLR DRCNT2

```

| | | | |
|------------|-----|---------|------|
| 11C9 3F 5F | 825 | CLR | REG3 |
| 11CB 81 | 826 | NODRP | RIS |
| | 827 | ; | |
| | 828 | ; | |
| | 829 | \$EJECT | |

```

830 ;
831 ;
832 ;
833 ; RPM CHECK
834 ;-----
835 ; This routine samples the hall effect line. At rates over 95 it stops
836 ; the motor every 3 magnets or one rev. At rates under 100 it stops
837 ; every magnet 1/3 rev.
838 RPM3 JMP RPM1
839 ;
840 ; Conditions
841 ;
842 RPMCK BRCIR 2,PSTAT1,NODRP ;DONT DO IF UNIT OFF
843 BRCIR 1,PSTAT1,NODRP ;DON'T DO IF MOTOR IS OFF
844 ;
845 ; Delay until circuit settles
846 ;
847 LDA MOTIM2 ;DON'T DO IF MOTOR NOT RUN ENOUGH
848 BNE RPM2 ;IF 2ND BYTE >0 THEN DO TEST
849 LDA MOTIM1 ;IF=0M THEN TEST 1ST BYTE FOR>$32
850 CMP #86 ;
851 BLO NODRP ;IF <THEN 86*.00816=.7017 SEC THEN DONT DO

```

```

852 ;
853 ; Test magnet
854 ;
855 RPM2 BRSET 0,PORTD,RPM3 . ;IF LINE HIGH THERE IS NO MAGNET
856 CLR ERRCNT
857 INC ITEST ;UPDATE MAG LOW COUNTER
858 LDA ITEST
859 CMP #255 ;IF MAG SENSED FOR 2 CONS. SEC
860 BEQ TIMERR ;THEN MAGNET OR ROTOR ERROR.
861 CMP #$01 ;IF ITEST > 1,
862 BHI NORPM ;SAME PULSE DETECTED SO DON'T UPDATE
863 ;
864 ;----- If this is model 324 add .376 to volume for every rev -----
865 ;
866 BRCIR 7,PSTAT1,OFMOT1 ;IF 224 DONT ADD
867 ;
868 JSR CALCV ;ELSE, CALCULATE VOLUME
869 ;
870 ;
871 ;
872 XYZV BRCIR 6,VOIM,OFMOT ;IF NOT IN VOLUME MODE OR
873 BRSET 0,VFLAG,OFMOT ;VOL TEST DISPLAY ENABLED

```

```

11FC CD 1634      874      JSR    VSEXP      ;GO UPDATE VOLUME NUMBERS
                   875      ;
                   876      ;
                   877      ;
                   878      ;
                   879      ;      Test 1/3 OR 1 REV  AND Turn off the motor
                   880      ;
                   881      OFMOT LDA DS1      ;IF DOSE!=0, CHECK DOSE
                   882      ORA DS2
11FF B6 55
1201 BA 56

1203 26 36      883      BNE    CMPDOS
                   884      ;
1205 3F 5B      885      OFMOT1 CLR REG1
1207 3C 61      886      INC REG5      ;CHECK IF THIS IS START UP
1209 B6 61      887      LDA REG5      ;IF SO IGNORE THE FIRST MAGNET
120B A1 02      888      CMP #02      ;AS IT COULD BE LESS THEN 1/3
120D 25 20      889      BLO NORPM      ;IT IS SO JUMP OVER SET CHECK
                   890      ;
120F 3A 61      891      DEC REG5      ;ITS NOT START SO KEEP FROM OVERFLOW
1211 3C 7A      892      INC THIRD1
1213 B6 7A      893      LDA THIRD1

```

```

1215 A1 03      894      CMP    #03
1217 26 0A      895      BNE    OFF90
1219 3F 7A      896      CLR    THIRD1
                        897      ;
                        898      ;      If drop flag not set then error
                        899      ;
121B 01 81 43   900      BRCLR 0,COUNT4,DROPER
121E 02 5F 40   901      BRSET 1,REG3,DROPER      ;EYES BLOCKED FOR 1 REV
                        902      ;
1221 3F 81      903      CLR    COUNT4
1223 B6 53      904      OFF90 LDA  DSPDG3      ;CHECK IF >95 RATE
1225 26 31      905      BNE    STP3      ;IF RATE >100 THEN STOP EVERY 3 SENSORS
1227 20 02      906      BRA    OFFM1
                        907      ;
1229 10 6D      908      OFFM BSET 0,PSTAT1      ;SET HOLD FLAG TO DISABLE MOTCK
122B 13 6D      909      OFFM1 BCLR 1,PSTAT1      ;MOTOR OFF PSTAT
122D 1B 01      910      BCLR 5,PORTB      ;MOTOR OFF
122F 81         911      NOREM RTS
                        912      ;
                        913      ;
                        914      ;
                        915      ;

```

```

916 ;
917 ;      Reset pulse status
918 ;
919 REM1 CIR      LTEST      ;CLEAR MAG LOW COUNTER
920 INC  ERCNT
921 LDA  ERCNT
922 CMP  #255      ;IF NO MAGNET IN 255X.008=2.04 SEC THEN ERR
923 BEQ  TIMERR
924 RTS
925 ;
926 ;
927 CMPDOS JSR    CPVTOD      ;CALL COMPARE VOL TO DOSE
928 CPX  #$FF      ;X IS FF IF VOL >= DOSE
929 BNE  OFMOIT1    ;IF X NOT FF, CONTINUE
930 ;              ;ELSE DOSERR
931 ;
932 DOSERR BSET 0,DOSER
933 LDA  PSTAT2
934 ORA  #$A0
935 AND  #$BD

```

5-7593

```

124A B7 6E      STA  PSTAT2
124C 20 DB      BRA  OFFM
936             ;
937             ;
938             ;
939             ;
940             TIMERR      BSET 7,ERRON      ;ENABLE ERROR FLAG
941             BSET 5,AIRM      ;ENABLE ALARM
1250 1A 6E      CLR  ERRCNT
942             JSR  OFFM
1252 3F 7F
943             RTS
1254 CD 1229
1257 81
944             ;
945             ;      Count every 3 magnet hits
946             ;
947             ;
948             STP3  INC  THIRDR      ;BUMP THIRD OF A REV COUNTER;
1258 3C 7B      LDA  THIRDR
949             CMP  #$03      ;IF THIRD ONE THEN STOP
125A B6 7B      BHS  OFFM1      ;DONT SET HOLD FLAG
125C A1 03
125E 24 CB      RTS
1260 81
950             ;
951             ;
952             ;      No drop error
953             ;
954             ;
955             ;
956             ;
957             DROPER      BSET 2,NDROP      ;SET FLOW ERROR FLAG

```



```

1263 1E 6E      BSET 7,ERRON      ;ENABLE DISPLAY ERROR MODE
1265 1A 6E      BSET 5,ALRM       ;ENABLE ALARM ERROR
1267 1D 6E      BCIR 6,VOLM      ;GET OUT OF VOLUME MODE
1269 11 71      BCIR 0,DFLAG     ;CLEAR DOSE MODE FLAG
126B 20 BC      BRA  OFFM
;
958             ;
964             ;               .125 ml/ 1/3 rev 8*.125=1
965             ;
126D 3C 58      CALCV INC DECML1  ;INCREMENT DECML COUNTER
126F B6 58      LDA  DECML1
1271 A1 08      CMP  #$08        ;IF DECML=8, TIME TO INCREMENT VS1, VS2
1273 26 1A      BNE  ENDCALC     ;ELSE, ENDCALC
;
1275 3F 58      CLR  DECML1     ;CLEAR DECML COUNTER
;
1277 A6 99      LDA  #$99       ;CHECK FOR FULL COUNTER
1279 B1 59      CMP  VS1
127B 26 04      BNE  CONCAL     ;IF VS1 != 99, CONTINUE CALCULATION
127D B1 5A      CMP  VS2       ;ELSE, CHECK VS2
127F 27 0E      BEQ  ENDCALC    ;IF VS1 & VS2 = 99, END CALCULATION
;
1281 B6 59      CONCAL LDA VS1  ;INCREMENT VS1

```

| | | | |
|------------|-----|-------------|------------------------|
| 1283 AB 01 | 980 | ADD #S01 | |
| 1285 8D | 981 | DAA | ;DECIMAL ADJUST RESULT |
| 1286 B7 59 | 982 | STA VS1 | ;STORE VS1 |
| | 983 | ; | |
| 1288 B6 5A | 984 | LDA VS2 | ;ADD VS2 |
| 128A A9 00 | 985 | ADC #S00 | |
| 128C 8D | 986 | DAA | ;DECIMAL ADJUST RESULT |
| 128D B7 5A | 987 | STA VS2 | ;STORE RESULT |
| 128F 81 | 988 | ENDCALC RTS | |
| | 989 | ; | |
| | 990 | \$EJECT | |

```

991 ;
992 ;          TIMING THE MOTOR
993 ;-----
994 ; This routine compares the calculated vs. timed motor on time and restarts
995 ; the motor every time they are equal.
996 ; The motor is turned off in the rpm routine>
997 ; the calculated values are generated by the UPDATE subroutine.
998 ;
999 ; Variables: MOTIM1=INCREMENTED EACH TIME
1000 ;           MOTIM2=INCREMENTED EACH TIME MOTIM1 OVERFLOWS
1001 ;           TIMLO= CALCULATED FIRST BYTE OF TIME ON GENERATED IN UPDATE
1002 ;           TIMHI= CALCULATED SECOND BYTE OF TIME
1003 ;
1004 ;           Conditions
1005 ;
1006 ;           MOTCK BRCLR 2, PSTAT1,END1 ;DONT DO IF UNIT OFF
1007 ;           BRSET 0, PSTAT1,END1 ;DON'T DO IF HOLD
1008 ;           BRSET 7, ERRON,END1 ;DON'T DO IF ERROR MODE
1009 ;           BRSET 3, PSTAT1,END1 ;DONT DO IF START MODE
1010 ;
1011 ;           ABOVE LINES DELETED TO ALLOW MOTOR TO
1012 ;           PUMP WITH LO BAT ERROR

```

```

1013 ;
1014 ;
1015 ; Check if time running equals time calculated
1016 ;
1017 CLC
1018 LDA MOTIM1
1019 ADD #$01
1020 STA MOTIM1
1021 BCS INCM2
1022 LDA TIMLO
1023 CMP MOTIM1
1024 BNE END1
1025 LDA TIMHI
1026 CMP MOTIM2
1027 BEQ MOTON
1028 RTS
1029 ;
1030 INCM2 CLC
1031 INC MOTIM2
1032 JMP BACK
1033 ;
1034 ; Times are equal so turn motor on again


```

3NSDOCID: <EP_____0327209A2_I_>

```

1045 ;
1046 ;
1047 ; SET CHECK
1048 ;-----
1049 ; Conditions
1050 ;
1051 SETCK BRCIR 2,PSTAT1,NOSCK ;DONT DO IF UNIT OFF
1052 BRSET 3,PSTAT1,NOSCK ;DONT DO IF UNIT IS IN START MODE
1053 BRSET 0,PSTAT1,NOSCK ;DONT DO IF UNIT IS IN HOLD
1054 ;
1055 ; CHECK FOR PROPER SET PLACEMENT
1056 ;
1057 BRCIR 3,PORTC,NOSCK ;IF SET PRESENT, RETURN
1058 ; ELSE SET ERROR
1059 BSET 4,NOSET ;ENABLE SET ERROR FLAG
1060 BSET 7,ERRON ;ENABLE ERROR MODE
1061 BSET 5,ALRM ;ENABLE AUDIO ALARM
1062 BCIR 6,VOLM ;GET OUT OF VOLUME MODE
1063 BCIR 0,DFLAG ;CLEAR DOSE MODE FLAG
1064 BCIR 1,PSTAT1 ;CLEAR MOTOR ON FLAG
1065 BCIR 5,PORTB ;TURN OFF MOTOR
1066 BSET 0,PSTAT1 ;SET HOLD FLAG TO DISABLE MOTCK

```

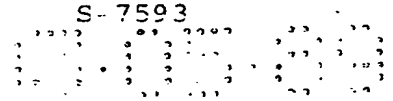

| | | | |
|-------|---|-----------|---------|
| 1067. | ; | NOSCK R1S | |
| 1068 | | | 12D9 81 |
| 1069 | | | |
| 1070 | | | |
| 1071 | | | |
| 1072 | | | |
| 1073 | | \$EJECT | |

ISDOCID: <EP_____0327209A2_I_>


```

12FB B7 82      1096 STA LSTRAT      ;STORE RATE POSITION FOR NEXT TIME
                  1097 ;
12FD A1 01      1098 CMP #01          ;BEGIN DECODING DEC POSITIONS
12FF 27 64      1099 BEQ DECR1
1301 A1 07      1100 CMP #07
1303 27 60      1101 BEQ DECR1
1305 A1 08      1102 CMP #08
1307 27 5C      1103 BEQ DECR1
1309 A1 0E      1104 CMP #0E
130B 27 58      1105 BEQ DECR1
                  1106 ;
130D A1 02      1107 CMP #02          ;BEGIN DECODING INC POSITIONS
130F 27 0D      1108 BEQ INCR1
1311 A1 04      1109 CMP #04
1313 27 09      1110 BEQ INCR1
1315 A1 0B      1111 CMP #0B
1317 27 05      1112 BEQ INCR1
1319 A1 0D      1113 CMP #0D
131B 27 01      1114 BEQ INCR1
131D 81          1115 OVER RTS
                  1116
                  1117 ;          test for increase by one

```

```

1118 ;
1119 INCR1 BC1R 5,ZEROST ;CLEAR DISPLAY ZERO FLAG
1120 BSET 3,PSTAT1 ;MAKE START ON
1121 BSET 4,PSTAT1 ;MAKE DISPLAY ON
1122 LDA #S03 ;IF 300, RETURN
1123 CMP DSPDG3
1124 BEQ OVER
1125 ;
1126 BRC1R 7,PSTAT1,UPFIVE ;IF MODEL 224 GO UP BY 5 ONLY
1127 ;
1128 ;IF RATE >= 50, INCREMENT BY 5
1129 LDA DSPDG3 ;
1130 BNE UPFIVE ;IF DIGIT 3 IS NOT 0, INCREMENT BY 5
1131 LDA DSPDG2
1132 CMP #S05
1133 BHS UPFIVE ;IF DIGIT 2 IS >= 5, INCREMENT BY 5
1134 ;
1135 UPONE CLC ;ELSE INCREMENT BY 1
1136 LDA #F7 ;TEST FOR 9
1137 ADD DSPDG1
1138 BCS DEC91 ;IF 9 MAKE 0 AND INC NEXT DIGIT
1139 INC DSPDG1 ;IF NOT THEN INC BY ONE

```

3NSDOCID: <EP_____0327209A2_I_>

```

1162 ;
1163 ; Increment digit 3 by one
1164 ;
1165 INC3 INC DSPDG3
1166 BCLR 7,ZEROST
1167 DEC3 JMP UPDATE
1168 ;
1169 ; Do a rate decrease
1170 ;
1171 ; test for decrement by 5
1172 ;
1173 DEC3 BSET 3,PSTAT1 ;MAKE START ON
1174 BSET 4,PSTAT1 ;MAKE DISPLAY ON
1175 BRSET 5,ZEROST,OVER ;IF DISPLAY IS 0 THEN DONT DECR
1176 LDA DSPDG3 ;PREPARE ACCA FOR ZERO TEST OF
1177 ORA DSPDG2 ;DSPDG2 AND DSPDG3
1178 BRSET 7,PSTAT1,CH324 ;IF 324 GOTO CH324, ELSE
1179 CH224 BNE DECR5 ;IF 224 AND DIGIT 2 & 3 ARE NOT ZERO, DECR5
1180 RTS ;ELSE DONT DECR
1181 ;
1182 CH324 BNE DECCN ;IF DIGIT2 & 3 ARE NOT ZERO, CONTINUE
1183 LDA DSPDG1 ;ELSE, TEST DIGIT 1 FOR 1

```

5-7593

```

137A A1 01      1184      CMP      #01
137C 26 13      1185      BNE      DOWN1      ;IF NOT 1,DECR BY 1
137E 81          1186      RTS
                  1187      ;
137F B6 53      1188      DECCN LDA  DSPDG3      ;IF DIGIT3 NOT ZERO,
1381 26 54      1189      BNE      DECR5      ;DECR BY 5
1383 B6 52      1190      LDA      DSPDG2      ;IF DIGIT2 > 5,
1385 A1 05      1191      CMP      #05      ;
1387 22 4E      1192      BHI      DECR5      ;DECR BY 5
1389 25 06      1193      BLO      DOWN1      ;ELSEIF DIGIT2 < 5 DECR BY 1
138B B6 51      1194      LDA      DSPDG1      ;ELSE TEST DIGIT1
138D A1 00      1195      CMP      #00      ;IF DIGIT1 NOT ZERO,
138F 26 46      1196      BNE      DECR5      ;DECR BY 5
                  1197      ;
                  1198      ;
                  1199      ;      DECREMENT BY 1
                  1200      ;
1391 B6 51      1201      DOWN1      LDA      DSPDG1      ;GET DIGIT 1
1393 27 09      1202      BEQ      TET2      ;IF ZERO THEN MAKE 9
1395 A1 01      1203      CMP      #01
1397 27 34      1204      BEQ      CHK300
1399 3A 51      1205      DEC2      DEC      DSPDG1      ;IF NOT THEN JUST DECREMENT

```

```

139B CC 1756      1206      JMP      UPDATE      ;GO UPDATE THE RATE
1207      ;
139E B6 52      1208      TET2 LDA DSPDG2
13A0 26 0B      1209      BNE INC9A
13A2 3A 53      1210      DEC DSPDG3
13A4 A6 09      1211      LDA # $09
13A6 B7 51      1212      STA DSPDG1
13A8 B7 52      1213      STA DSPDG2
13AA CC 1756      1214      JMP      UPDATE
1215      ;
13AD 3A 52      1216      INC9A DEC DSPDG2      ;IF NOT DECREMENT BY ONE
13AF A6 09      1217      LDA # $09      ;MAKE FIRST DIGIT 9
13B1 B7 51      1218      STA DSPDG1
13B3 CC 1756      1219      JMP      UPDATE      ;GO UPDATE THE DISPLAY
1220      ;
13B6 A6 03      1221      GO300 LDA # $03      ;CHANGE DIGIT 3 TO 3
13B8 B7 53      1222      STA DSPDG3
13BA A6 00      1223      LDA # $00      ;MAKE DISPLAY 300
13BC B7 52      1224      STA DSPDG2
13BE B7 51      1225      STA DSPDG1
1226      ;      CLR ZEROST
1227      LDA SPEED1

```

| | | | | | | | |
|--------------|------|-----------|--------|--------|--|--|-----------------------------------|
| 13C2 A1 40 | 1228 | CMP | #\$40 | | | | |
| 13C4 25 04 | 1229 | BLO | SPDY1 | | | | |
| 13C6 A6 41 | 1230 | LDA | #\$41 | | | | |
| 13C8 B7 78 | 1231 | STA | SPEED1 | | | | |
| 13CA CC 1756 | 1232 | SPDY1 JMP | UPDATE | | | | |
| | 1233 | ; | | | | | |
| 13CD B6 52 | 1234 | CHK300 | LDA | DSPDG2 | | | ;IF DIGIT 2 IS ZERO MAKE 300 |
| 13CF 26 C8 | 1235 | BNE | DEC2 | | | | |
| 13D1 B6 53 | 1236 | LDA | DSPDG3 | | | | |
| 13D3 26 C4 | 1237 | BNE | DEC2 | | | | |
| 13D5 20 DF | 1238 | BRA | GO300 | | | | |
| | 1239 | ; | | | | | |
| | 1240 | ; | | | | | change rate by 5 |
| | 1241 | ; | | | | | |
| 13D7 4F | 1242 | DECR5 | CLRA | | | | |
| 13D8 B1 53 | 1243 | CMP | DSPDG3 | | | | ;IF DIGIT 2,3 ARE 0 THEN MAKE 000 |
| 13DA 26 0D | 1244 | BNE | TOG2 | | | | |
| 13DC B1 52 | 1245 | CMP | DSPDG2 | | | | |
| 13DE 26 09 | 1246 | BNE | TOG2 | | | | |
| 13E0 B7 51 | 1247 | STA | DSPDG1 | | | | |
| 13E2 A6 03 | 1248 | LDA | #\$03 | | | | |
| 13E4 B7 53 | 1249 | STA | DSPDG3 | | | | |

| | | | | | |
|--------------|------|---|----------|--------|----------------------|
| 13E6 CC 1756 | 1250 | ; | CLR | ZEROST | |
| 13E9 CD 174B | 1251 | | JMP | UPDATE | |
| 13EC A6 00 | 1252 | | TOG2 JSR | TOGGLE | ;GO MAKE 0 OR 5 |
| 13EE B1 51 | 1253 | | LDA | #\$00 | |
| 13F0 27 19 | 1254 | | CMP | DSPDG1 | |
| | 1255 | | BEQ | INCD3 | |
| | 1256 | ; | | | |
| | 1257 | ; | | | |
| | 1258 | ; | | | |
| 13F2 B6 52 | 1259 | | LDA | DSPDG2 | ;GET POINTER |
| 13F4 27 05 | 1260 | | BEQ | INC9 | ;IF 0 THEN MAKE 9 |
| 13F6 3A 52 | 1261 | | DEC | DSPDG2 | ;IF>0 THEN SUB 1 |
| 13F8 CC 1756 | 1262 | | SETZ | JMP | UPDATE |
| 13FB B6 53 | 1263 | | INC9 | LDA | DSPDG3 |
| 13FD 27 0C | 1264 | | BEQ | INCD3 | |
| 13FF A6 09 | 1265 | | LDA | #\$09 | ;RESET TO 9 |
| 1401 B7 52 | 1266 | | STA | DSPDG2 | |
| | 1267 | ; | | | |
| | 1268 | ; | | | |
| | 1269 | ; | | | |
| | 1270 | | DECD3 | LDA | DSPDG3 |
| 1403 B6 53 | 1271 | | BEQ | INCD3 | ;GET DIGIT 3 POINTER |
| 1405 27 04 | | | | | ;IF 0 THEN MAKE 2 |

EP 0 327 209 A2

S-1
75593

1407 3A 53
1409 1F 74
140B CC 1756

| | | |
|------|---------|------------|
| 1272 | DEC | DSPOG3 |
| 1273 | BCLR | 7,ZEROST |
| 1274 | INCD3 | JMP UPDATE |
| 1275 | : | |
| 1276 | : | |
| 1277 | \$EJECT | |

```
!>0 SO SUB 1
```

EP 0 327 209 A2

S-7593

| | | |
|--------------|-------|--|
| 1278 | ; | |
| 1279 | ; | |
| 1280 | ; | OFF CHECK |
| 1281 | ; | |
| 1282 | ; | Conditions |
| 1283 | ; | |
| 1284 | OFFCK | BRCLR 2, PSTAT1, TSTCK ; IF OFF, CHECK IF TIME FOR TEST MODE |
| 1285 | BRCLR | 2, PORTC, GOFF ; IF OFF BUTTON PRESSED GO OFF |
| 1414 81 | NOOFF | RTS |
| 1286 | ; | |
| 1287 | ; | |
| 1288 | ; | Do off |
| 1289 | ; | |
| 1290 | GOFF | BRSET 0, LOWBAT, KIPLP1 ; IF LOW BATT KILL POWER |
| 1418 3F 6E | CLR | PSTAT2 ; NO ALARMS OR OPTIONAL MODES |
| 141A 3F 71 | CLR | DELAG ; CLEAR DOSE FLAG |
| 141C 3F 80 | CLR | DSP1ST ; CLEAR DISPLAY TEST FLAG |
| 1293 | | |
| 1294 | ; | |
| 141E A6 19 | LDA | #\$19 ; SET PSTAT1 FOR NEXT TURN-ON |
| 1295 | | |
| 1420 B7 6D | STA | PSTAT1 |
| 1296 | | |
| 1422 A6 1F | LDA | #\$1F ; KILL DISPLAY AND MOTOR 0001 1111 OR |
| 1297 | | |
| 1424 B7 01 | STA | PORTB |
| 1298 | | |
| 1426 CD 1948 | JSR | TM1ST ; CLEAR ALL TIMERS FOR 24HR KEEPING |
| 1299 | | |

```

1429 3F 45      1300      CLR      T25
142B 3F 55      1301      CLR      DS1      ;CLEAR DS1 AND DS2 FOR USE AS
142D 3F 56      1302      CLR      DS2      ;3 SEC TEST TIMER.
142F 8F          1303      LOOFP      WAIT      ; GO SLEEP
1430 20 FD      1304      BRA      LOOFP
1305      ;
1306      ;      wake up by on button (external int)
1307      ;
1432 CD 1106     1308      JSR      TIMEON      ;RESET COUNTER TIMER
1435 81          1309      RTS
1310      ;
1436 9B          1311      KILLP1      SEI      ;SET INT MASK TO PREVENT TURN-ON
1437 A6 0F      1312      LDA      #$0F      ;DISABLE PUMP
1439 B7 0B      1313      STA      PORTB
143B 20 F9      1314      BRA      KILLP1      ;LOOP UNTILL POWER DISSIPATES
1315      ;
1316      ;      TEST CHECK
1317      ;
143D 04 02 03    1318      TSTCK BRSET 2,PORTC,TSTOUT ;IF OFF BUTTON NOT PUSHED, RETURN
1440 05 03 05    1319      BRCLR 2,PORTD,TSTCNT ;ELSE, IF HOLD BUTTON PRESSED INCREMENT
1320      ;      ; YES COUNTER
1443 3F 55      1321      TSTOUT CLR DS1

```

EP 0 327 209 A2

S-7593

| | | | |
|--------------|------|-----------------|---------------------------------|
| 1445 3F 56 | 1322 | CLR DS2 | |
| 1447 81 | 1323 | RTS | |
| | 1324 | ; | |
| 1448 3C 55 | 1325 | TSTCNT INC DS1 | : INCREMENT TEST TIMER |
| 144A 26 02 | 1326 | BNE TST1 | |
| 144C 3C 56 | 1327 | INC DS2 | |
| 144E B6 56 | 1328 | TST1 LDA DS2 | : IF 3 SEC ELAPSED, GOTO PMPTST |
| 1450 A1 02 | 1329 | CMP #502 | |
| 1452 24 01 | 1330 | BHS PTST | |
| 1454 81 | 1331 | RTS | : ELSE RETURN |
| 1455 CC 1C00 | 1332 | PTST JMP PMPTST | |
| | 1333 | ; | |
| | 1334 | \$EJECT | |

| | | | |
|---------------|------|---|--|
| 1458 05 6D 30 | 1335 | : | |
| 145B 01 80 2D | 1336 | : | |
| 145E 0F 6D 2D | 1337 | : | |
| 1461 3C 79 | 1338 | : | |
| 1463 B6 79 | 1339 | : | |
| 1465 A1 78 | 1340 | : | |
| 1467 22 4F | 1341 | : | |
| 1469 A1 64 | 1342 | : | |
| 146B 22 7B | 1343 | : | |
| 146D A1 50 | 1344 | : | |
| 146F 22 66 | 1345 | : | |
| 1471 A1 3C | 1346 | : | |
| 1473 22 5E | 1347 | : | |
| 1475 A1 28 | 1348 | : | |
| 1477 22 49 | 1349 | : | |
| 1479 A1 14 | 1350 | : | |
| 147B 22 22 | 1351 | : | |
| | 1352 | : | |
| | 1353 | : | |
| | 1354 | : | |
| | 1355 | : | |
| | 1356 | : | |

| | | TEST MODE |
|-------------------------------|--------------------------------|-----------|
| TESTED BRCLR 2, PSTAT1, NOTST | : DONT DO IF UNIT OFF | |
| BRCLR 0, DSPST, NOTST | : DONT DO IF TEST FLAG NOT SET | |
| BRCLR 7, PSTAT1, TSTD2 | : IF 224, DO SHORT TEST MODE | |
| INC VOLFTM | : DO FOR .033*60=1.98 SEC | |
| LDA VOLFTM | | |
| CMP #120 | : GO CLEAR TEST MODE | |
| BHI CLTEST | | |
| CMP #100 | : GO DISPLAY DOSE NUMBER | |
| BHI DOSDIS | | |
| CMP #80 | : GO DISPLAY DOSE | |
| BHI DOSTST | | |
| CMP #60 | : GO DISPLAY VOLUME NUMBER | |
| BHI VOLDS | | |
| CMP #40 | : GO DISPLAY VOL | |
| BHI VOLNST | | |
| CMP #20 | : GO KILL ALARM FIRST | |
| BHI CLATST | | |

EP 0 327 209 A2

```

147D 1E 01      1357      CONTST      BSET 7,PORTB      ;ALARM ON
147F 1C 01      1358      BSET 6,PORTB      ;CLEAR LED ON
1481 1C 6E      1359      BSET 6,VOIM      ;SET VOIM FLAG TO ENABL 4 DIGITS
1483 3F 4D      1360      CIR  SAMPL1      ;MAKE DISPLAY "8888"
1485 3F 4E      1361      CIR  SAMPL2
1487 3F 4F      1362      CIR  SAMPL3
1489 3F 50      1363      CIR  SAMPL4
148B 1D 03      1364      NOTST BCIR 6,PORTD      ;TURN IR SOURCE OFF
148D 81         1365      RTS
1366           1366      ;
1367           1367      ; SHORT TEST FOR K224
1368           1368      ;
148E 1C 6E      1369      TSTD2 BSET 6,PSTAT2      ;VOLUME FLAG SET TO ENABLE 4 DIGITS
1490 3C 79      1370      INC  VOLTIM      ;DO FOR .033*60=1.98 SEC
1492 B6 79      1371      LDA  VOLTIM
1494 A1 3C      1372      CMP  #$3C
1496 22 20      1373      BHI  CLTEST
1498 A1 14      1374      CMP  #$14
149A 22 03      1375      BHI  CLATST
149C CC 147D    1376      JMP  CONTST
1377           1377      ;
1378           1378

```

```

149F 1F 01      1379      CLATST      BC1R 7,PORTB      ;TURN ALARM OFF
14A1 1D 01      1380      BC1R 6,PORTB      ;LED OFF
1381      ;
1382      ;      test sensor eyes
1383      ;
1384      BSET 6,PORTD      ;TURN IR SOURCE ON
14A3 1C 03      1385      NOP      ;ADD 8 USEC DELAY TO ALLOW IR
14A5 9D      1386      NOP      ;SENSOR TO SETTLE
14A6 9D      1387      NOP
14A7 9D      1388      NOP
14A8 9D      1389      BRSET 5,PORTC,NOTST      ;IF RECIEVER LOW EYES BLOCKED THEN ERROR
14A9 0A 02 DF      1390      BC1R 6,PORTD      ;TURN OFF IR SOURCE
14AC 1D 03      1391      BSET 5,PSTAT2      ;ENABLE ALARM
14AE 1A 6E      1392      BSET 2,NDROP      ;ENABLE FLOW ERROR FLAG
14B0 14 6E      1393      BSET 7,PSTAT2      ;ENABLE ERROR FLAG
14B2 1E 6E      1394      BC1R 7,PORTB      ;TURN ALARM OFF
14B4 1F 01      1395      BC1R 3,PSTAT1      ;CLEAR START FLAG
14B6 17 6D      1396      ;
14B8 CD 1756      1397      CLITEST      JSR UPDATE      ;RESTORE DISPLAY
14BB 3F 80      1398      CLR DSPTST      ;CLEAR TEST FLAG
14BD 1D 6E      1399      BC1R 6,VOLM      ;CLEAR VOLUME FLAG FOR 224
14BF 3F 79      1400      CLR VOLTIM      ;CLEAR TIMER

```

5 7 5 9 3

```

14C1 81      1401      RTS
                ;
14C2 A6 F8    1402      VOLTST      LDA  #$F8      ;
14C4 B7 50    1404      STA  SAMPLE4
14C6 A6 9C    1405      LDA  #$9C
14C8 B7 4F    1406      STA  SAMPLE3
14CA A6 C8    1407      LDA  #$C8
14CC B7 4E    1408      STA  SAMPLE2
14CE A6 FF    1409      LDA  #$FF
14D0 B7 4D    1410      STA  SAMPLE1
14D2 81      1411      RTS
                ;
14D3 CD 1634  1413      VOLDS JSR  VSEXP
14D6 81      1414      RTS
                ;
14D7 A6 B0    1416      DOSIST      LDA  #$B0
14D9 B7 50    1417      STA  SAMPLE4
14DB A6 92    1418      LDA  #$92
14DD B7 4F    1419      STA  SAMPLE3
14DF A6 9C    1420      LDA  #$9C
14E1 B7 4E    1421      STA  SAMPLE2
14E3 A6 8C    1422      LDA  #$8C
    
```


| | | | | | | | | | |
|--------------|------|--------|--------|-------|--|--|--|--|--------------------------|
| 14E5 B7 4D | 1423 | STA | SAMPL1 | | | | | | |
| 14E7 81 | 1424 | RTS | | | | | | | |
| | 1425 | | | | | | | | |
| | | | | | | | | | |
| 14E8 CD 1745 | 1426 | DOSDIS | JSR | DSEXP | | | | | ;DISPLAY PROGRAMMED DOSE |
| 14EB 81 | 1427 | RTS | | | | | | | |
| | 1428 | | | | | | | | |
| | 1429 | | | | | | | | |
| | 1430 | | | | | | | | |
| | | | | | | | | | |
| | 1431 | | | | | | | | \$EJECT |

```

1432 ;
1433 ; HOLD CHECK
1434 ;-----
1435 ; Conditions
1436 ;
1437 HOLDCK BRCLR 2,PSTAT1,NOHOLD ;IF OFF DONT DO
1438 BRSET 0,DSPTST,NOHOLD ;DONT DO IF TEST MODE
1439 BRSET 0,CWF,NOHOLD ;DONT DO IF INC OR DEC
1440 BRSET 0,CCWF,NOHOLD ;BUTTON IS PRESSED
1441 ;
1442 BRSET 2,PORID,NOH1 ;IS HOLD BUTTON DOWN RETURN IF NOT
1443 ;
1444 ; Hold button down
1445 ;
1446 ; Debounce switch
1447 ;
1448 BRSET 0,HTEST,NOHOLD ;WE MUST BE SEEING SAME PULSE SO IGNORE
1449 BSET 0,HTEST ;IF NOT SAME SET THE FLAG
1450 ;
1451 ; CLEAR HOLD TIME COUNTERS
1452 ;
1453 CLR HLDIM1
1500 3F 4B

```

| | | | | |
|---------------|------|--------------------------|--|--|
| 1502 3F 4C | 1454 | CIR | HIDTM2 | |
| | 1455 | | | ; |
| | 1456 | | Check for alarm conditions | ; |
| | 1457 | | | ; |
| 1504 00 6E 2F | 1458 | BRSET 0, LOWBAT, HLD BAT | | ; IF LOW BAT, STOP ALARM AND MOTOR. |
| 1507 00 70 37 | 1459 | BRSET 0, DOSER, STDSE | | ; IF DOSE DEL FLAG SET, STOP DOSE MESS. |
| 150A 0A 6E 2F | 1460 | BRSET 5, ALPM, NOAUDO | | ; IF WE ARE ALARMING THEN KILL AUDIO |
| 150D 0E 6E 3F | 1461 | BRSET 7, ERRON, NOMESS | | ; IF WE ARE FLASHING ERROR MESSAGE THEN GO RUN |
| 1510 0C 6E 34 | 1462 | BRSET 6, VOIM, NOVOIM | | ; IF IN VOL GO RESTORE RATE |
| 1513 00 71 35 | 1463 | BRSET 0, DFLAG, NODOS | | ; IF DOSE MODE, GO RESTORE RATE |
| | 1464 | | | ; |
| | 1465 | | Are we in hold now run if yes, hold if not | |
| 1516 06 6D 3E | 1466 | BRSET 3, PSTAT1, RUN | | ; IF IN START MODE THEN DO RUN |
| | 1467 | | | ; |
| 1519 00 6D 3B | 1468 | BRSET 0, PSTAT1, RUN | | ; IF YES GO RUN |
| 151C 10 6D | 1469 | HOLD BSET 0, PSTAT1 | | ; IF NOT SET HOLD FLAG |
| 151E 1B 01 | 1470 | BCIR 5, PORTB | | ; MOTOR OFF |
| 1520 13 6D | 1471 | BCIR 1, PSTAT1 | | ; CLEAR MOTOR ON FLAG |
| | 1472 | | | ; |
| 1522 B6 02 | 1473 | LDA PORTC | | ; RE-INITIALIZE RATE KNOB |
| 1524 A4 03 | 1474 | AND #03 | | |
| 1526 B7 82 | 1475 | STA LSTRAT | | |

```

1476 ;
1477 CLR PSTAT2
1478 CLR DOSER ;CLR DOSE ERROR FLAG
1479 JSR UPDATE ;RESTORE DIGITS TO KILL SCROLLING DP
1480 JSR TIMRST
1481 NOHOLD RTS
1482 ;
1483 ; Hold button not down so reset the flag and return
1484 ;
1485 NOH1 BCLR 0,HTEST
1486 RTS
1487 ;
1488 ; LOW BAT HOLD
1489 ;
1490 ; kill motor
1491 ;
1492 HLDSTAT BCLR 1,PSTAT1 ;CLEAR MOTOR ON FLAG
1493 BCLR 5,PORTB ;TURN MOTOR OFF
1494 BSET 0,PSTAT1 ;SET HOLD FLAG
1495 ;
1496 ; AND
1497 ;

```

```

1498 ; Kill the audio
1499 ;
1500 NOAUDIO BCIR 5,AIRM ;RESET BEEPER FLAG
1501 BCIR 7,PORTB ;STOP AUDIO NOW
1502 RTS
1503 ;
1504 ; STOP DOSE DELIVIERED MESSAGE
1505 ;
1506 STDOSE BCIR 6,VOIM ;CLEAR ALARM FLAG
1507 BCIR 7,PORTB ;STOP ALARM
1508 BRA HOLD ;PUT PUMP IN HOLD
1509 ;
1510 ; Recovery from error KILL THE MESSAGE
1511 ;
1512 NOVOIM BCIR 6,VOIM
1513 BRA HOLD
1514 ;
1515 NODOS BCIR 0,DFLAG
1516 BRA HOLD
1517 ;
1518 ;
1519 NOMESS CLR PSTAT2

153C 1B 6E
153E 1F 01
1540 81

1541 1D 6E
1543 1F 01
1545 20 D5

1547 1D 6E
1549 20 D1

154B 11 71
154D 20 CD

154F 3F 6E

```

```

1551 CD 1756      JSR  UPDATE      ;RESTORE DISPLAY TO RATE
1554 CC 1557      JMP  RUN          ;GO RUN MOTOR
1522             ;
1523             ;
1524             ;
1525             ;----- RUN MOTOR PUMPING -----
1526             ;
1527             ;
1528 RUN          BCIR 3,PSTAT1      ;SET UNIT START OFF
1529             BSET 4,PSTAT1      ;SET DISPLAY ON
1530             ;
1531             LDA  DS1            ;CHECK IF DOSE DEL BEFORE TURNING
1532             ORA  DS2            ;MOTOR ON
1533             BEQ  RUN1           ;IF DOSE = 0, NO NEED TO CHECK
1534             JSR  CPVTOD        ;ELSE, CALL CMP VOL TO DOSE ROUTINE
1535             CPX  #$FF          ;X IS FF IF VOL >= DOSE
1536             BNE  RUN1         ;IF X != FF, OK TO START MOTOR
1537             ;
1538             BSET 0,DOSER        ;ELSE, SET DOSE DEL MESSAGE
1539             BSET 0,PSTAT1      ;SET HOLD FLAG TO DISABLE MOTCK
1540             BSET 7,ERRON       ;ENABLE ERROR FLAG
1541             BSET 5,ALRM        ;ENABLE ALARM

```

| | | | |
|------------|------|-----------------|-----------------------|
| 1570 1D 6E | 1542 | BCLR 6, VOIM | ;CLEAR VOL MODE FLAG |
| 1572 11 71 | 1543 | BCLR 0, DFLAG | ;CLEAR DOSE MODE FLAG |
| 1574 19 6E | 1544 | BCLR 4, NOSET | ;CLR NO SET FLAG |
| 1576 81 | 1545 | RTS | |
| | 1546 | | ; |
| | 1547 | | ; |
| 1577 3F 7D | 1548 | RUN1 CLR MOTIM1 | |
| 1579 3F 7E | 1549 | CLR MOTIM2 | |
| 157B 3F 81 | 1550 | CLR COUNT4 | |
| 157D 3F 63 | 1551 | CLR COUNT2 | |
| 157F 3F 5B | 1552 | CLR REG1 | |
| 1581 3F 5D | 1553 | CLR DRCNT1 | |
| 1583 3F 5E | 1554 | CLR DRCNT2 | |
| 1585 3F 5F | 1555 | CLR REG3 | |
| 1587 3F 60 | 1556 | CLR REG4 | |
| 1589 3F 61 | 1557 | CLR REG5 | |
| 158B 3F 7A | 1558 | CLR THIRDR1 | |
| 158D 3F 7B | 1559 | CLR THIRDR | |
| 158F 3F 7F | 1560 | CLR ERRCNT | |
| 1591 3F 6E | 1561 | CLR PSTAT2 | |
| | 1562 | LDA #\$10 | ; |
| | 1563 | STA PSTAT2 | ; |

```
1564 ;
1565 LDA DSPDG1 ;CHECK RATE VALUE
1566 ORA DSPDG2
1567 ORA DSPDG3
1568 ORA DSPDG4
1569 BNE RUNMOT IF RATE!=0, RUNMOT
1570 ; BSET 7,ERRON ;IT =000 SO SET ERROR BIT
1571 BSET 5,ALRM ;SET ALARM BIT
1572 BSET 0,PSTAT1
1573 ; BCLR 0,PSTAT1
1574 RTS ;AND RETURN
1575 RUNMOT BSET 1,PSTAT1 ;MOTOR ON FLAG
1576 BCLR 0,PSTAT1 ;RUN FLAG ON
1577 BSET 5,POR1B ;MOTOR ON
1578 RTS
1579 ;
1580
1581
1582 $EJECT
```


7593
S
7593

```

1583 ; DOSE CHECK
1584 ; -----
1585 ; Conditions
1586 ;
1587 DOSECK BRCIR 2,PSTAT1,NDOSE ;DON'T DO IF UNIT IS OFF
1588 BRSET 0,DSPTST,NDOSE ;DON'T DO IF TEST MODE
1589 BRCIR 7,PSTAT1,NDOSE ;DON'T DO IF UNIT IS 224
1590 BRSET 7,PSTAT2,NDOSE ;DON'T DO IF ERROR MODE ON
1591 ;
1592 BRSET 0,DFLAG,DOS1 ;IF DOSE FLAG SET THEN DOS1
1593 ;
1594 BRSET 0,CWF,NDOSE ;IF INC OR DEC FLAG PRESSED
1595 BRSET 0,CCWF,NDOSE ;DONT SET UP DOSE FEATURE
1596 ;
1597 BRCIR 1,PORTD,SETDOS ;IF DOSE BUTTON PUSHED THEN
1598 ;INITIALIZE DOSE MODE
1599 RTS ;ELSE RETURN
1600 ;
1601 DOS1 INC DOSECT
1602 LDA DOSECT
1603 CMP #30
1604 BHI DOS2

```

```

15CA CD 14D7      JSR  DOSTST
15CD 81           NDOSE RTS
1605              ;
1606              SETDOS      BSET 0,DFLAG      ;SET DOSE MODE FLAG
1607              BCIR 6,VOLM      ;CLEAR VOLUME FLAG
1608              ;
1609              BRCLR 0,PSTAT1,CONDO      ;IF PUMP IN HOLD,
1610              BSET 3,PSTAT1      ;SET START FLAG TO DISABLE FLASH
1611              BSET 4,PSTAT1      ;MAKE DISPLAY ON
1612              ;ELSE, JUST CONTINUE WITH DOSE
1613              CONDO CLR  DOSECT      ;CLEAR DOSE COUNTER
1614              BRA  DOS1      ;CONTINUE DOSE HANDLING
1615              ;
1616              DOS2 JSR  DSEXP
1617              LDA  DOSECT
1618              CMP  #120
1619              BLO  NDOSE
1620              JSR  UPDATE
1621              BCIR 0,DFLAG      ;NO DOSE MODE
1622              CLR  DOSECT
1623              RTS
1624              ;
1625              ;
1626              ;

```

\$EJECT

1627

```

1628 ;
1629 ; VOLUME DISPLAY CHECK
1630 ;-----
1631 ; Conditions
1632 ;
1633 VOLDIS BRCIR 2,PSTAT1,NVL ;DONT DO IF OFF
1634 BRSET 0,DSPTST,NVL ;DONT DO IF IN TEST MODE
1635 BRCIR 7,PSTAT1,NVL ;DONT DO IF MODEL 224
1636 ; BRSET 0,DFLAG,NVL ;DON'T DO IF IN DOSE MODE
1637 BRSET 7,ERRON,NVL ;DON'T DO IF ERROR MODE
1638 BRSET 0,CWF,NVL ;DONT DO IF INC OR DEC BUTTON
1639 BRSET 0,CCWF,NVL ;IS PRESSED
1640 ;
1641 LDA VS1 ;UPDATE CLEAR VOLUME LED.
1642 ORA VS2
1643 BEQ NOVLED
1644 BCIR 6,PORTB ;TURN VOLUME LED OFF IF VOLUME NOT ZERO
1645 BRA VOLBUT1
1646 NOVLED BSET 6,PORTB ;IF NO VOLUME TURN ON LED
1647 ;
1648 VOLBUT1 BRSET 6,VOLM,VOL1 ;IF VOL FLAG SET UPDATE DISPLAY
1649 BRCIR 3,PORTD,SETVOL ;IF VOL BUTTON PUSHED, INITIALIZE DISPLAY

```

```

1612 81      RTS
1650
1651      ;
1652      SETVOL      BSET 6,VOLM      ;SET VOLUME DISPLAY MODE FLAG
1653      BCIR 0,DFLAG      ;CLEAR DOSE FLAG
1654      CIR  VOLTIM      ;CLEAR VOLUME TIMER
1655      BRCIR 0,PSTAT1,VOL1      ;IF PUMP IN HOLD MODE,
1656      BSET 3,PSTAT1      ;SET START FLAG TO DISABLE FLASH
1657      BSET 4,PSTAT1      ;MAKE DISPLAY ON
1658      ;ELSE, JUST CONTINUE VOLUME
1659      ;
1660      VOL1 INC  VOLTIM      ;INCREMENT VOLTIM
1661      LDA  VOLTIM
1662      CMP  #30      ;IF ELAPSED TIME IS
1663      BHI  VOL2      ; > 1.00 S, DISPLAY VOLUME
1664      JSR  VOLTST      ;ELSE DISPLAY VOLTST
1665      BSET 0,VFLAG      ;SET VOL TEST DISPLAY FLAG
1666      NVL  RTS
1667      ;
1668      VOL2 CMP  #120      ;CHECK ELAPSED TIME, ACCA LOADED
1669      BHI  ENDVOL      ;WITH VOLTIM IN VOL1 ROUTINE
1670      ;IF TIME > 5 SEC, END VOL MODE
1671      ;

```

```
1632 11 75      BC1R 0,VFLAG      ;CLEAR VOL TEST DISPLAY FLAG
1672
1673      ;
1674      VSEXP LDX #VS1      ;SET POINTER TO VS1 ADDRESS
1675      JSR BCDEXP
1676      RTS
1677      ;
1678      ENDVOL      JSR UPDATE      ;RETURN DISPLAY TO RATE
1679      BC1R 6,VOLM
1680      RTS
1681      ;
1682      ;
1683      ;
1684      ;
1685      $EJECT
```

```

1686 ;
1687 ;
1688 ;
1689 ; CHECK CLEAR VOLUME BUTTON
1690 ;-----
1691 ; Conditions
1692 ;
1693 CKCLR BRCLR 2,PSTAT1,NOCLR ;DONT DO IF OFF
1694 BRSET 0,DSPTST,NOCLR ;DONT DO IF IN TEST MODE
1695 BRCLR 7,PSTAT1,NOCLR ;DONT DO IF MODEL 224
1696 BRSET 0,CKBEP,CNALRM ;IF ALARM ALREADY SET, HANDLE BEEP
1697 ;
1698 ; Has volume clr button been pressed
1699 ;
1700 CLR BRSET 4,PORTD,NC1R1 ;IF CLR BUTTON NOT PRESSED RETURN
1701 ;
1702 ; DEBOUNCE SWITCH
1703 BRSET 1,HTEST,NOCLR ;IF BUTTON SENSED LAST CYCLE, MUST BE
1704 ; SAME PULSE SO IGNORE
1705 BSET 1,HTEST ;ELSE SET FLAG FOR NEXT CYCLE
1706 ;
1707 BRCLR 6,VOLM,NOCLR ;IF NOT VOLUME MODE, RETURN

```

```

1708 ;
1709 LDA VOLTIM ;CHECK VOLTIM
1710 CMP #31 ;IF 'VOL' DISPLAY STILL ON
1711 BLS NOCLR ;DON'T CLEAR YET
1712 ;
1713 CHECK IF OK TO CLEAR VOLUME
1714 ;
1715 LDA DS1 ;IF DOSE = 0, CONTINUE CLEAR
1716 ORA DS2
1717 BEQ CLR1
1718 ;
1719 JSR CPVTOD ;IF DOSE != 0 AND DOSE DEL.,
1720 CPX #$FF
1721 BEQ CLR1 ;IF X = FF, CONTINUE WITH CLEAR VOL
1722 ; ;ELSE, ENABLE SHORT ALARM
1723 CKALRM BSET 0,CKBEP ;SET CHECK ALARM BIT TO ENABLE BEEP
1724 CLR CKTIM ;INITIALIZE BEEP TIME
1725 BSET 7,PORTB ;TURN ON ALARM
1726 ;
1727 CVALRM INC CKTIM ;INCREMENT BEEP TIMER
1728 LDA CKTIM ;CHECK IF TIME TO STOP BEEP
1729 CMP #$05

```


| | | | |
|--------------|------|--------------------|-------------------------------|
| 1676 25 04 | 1730 | BLO CONBP | ;IF CKTIM < AA, CONTINUE BEEP |
| 1678 1F 01 | 1731 | BCLR 7,PORTB | ;ELSE, TURN OFF BEEP |
| 167A 11 83 | 1732 | BCLR 0,CKBEP | ;CLEAR CK BEEP FLAG |
| 167C 81 | 1733 | CONBP RTS | |
| | 1734 | ; | |
| | 1735 | ; CLEAR VOLUME | |
| | 1736 | ; | |
| 167D CD 1941 | 1737 | CIR1 JSR CLVOL | ;CLEAR VOLUME |
| 1680 CC 1634 | 1738 | JMP VSEXP | |
| | 1739 | ; | |
| 1683 81 | 1740 | NOCIR RTS | |
| | 1741 | ; | |
| 1684 13 73 | 1742 | NCIR1 BCLR 1,HTEST | |
| 1686 81 | 1743 | RTS | |
| | 1744 | ; | |
| | 1745 | \$EJECT | |

```

1746 ;
1747
1748 ;
1749 ;
1750 ;      Inputting the rate (model 324)
1751 ;      -----
1752 ;
1753 ;      Conditions
1754 ;
1687 05 6D 23 INRATE2  BRCLR 2,PSTAT1,OVER1  ;IF OFF DONT DO
168A 00 6E 20  BRSET 0,LOWBAT,OVER1  ;DONT DO IF LOWBAT
168D 0E 6E 1D  BRSET 7,ERRON,OVER1  ;DONT DO IF ERROR MODE
1690 0F 6D 1A  BRCLR 7,PSTAT1,OVER1  ;IF MODEL 224 DONT DO
1693 00 80 17  BRSET 0,DSPIST,OVER1  ;DONT DO IF IN TEST MODE
1696 0C 6E 14  BRSET 6,VOLM,OVER1  ;DONT DO IF IN VOLUME DISPLAY
1699 06 6D 03  BRSET 3,PSTAT1,CRIE  ;IF IN START MODE DO NOW
169C 01 6D 0E  BRCLR 0,PSTAT1,OVER1  ;DON'T DO IN RUN MODE
1763 ;
1764 ;      Check for rate change
1765 ;
169F 03 02 12  CRIE  BRCLR 1,PORTC,CCW1  ;IF BUTTON DOWN DO DECREASE
16A2 3F 77  CLR  CCWF  ;CLEAR FLAG IF NOT DOWN

```

```

16A4 01 02 07      1768      BRCLR 0,PORTC,CW1 ;IF BUTTON DOWN GO INCREASE
16A7 3F 76         1769      CLR  CWF      ;CLEAR FLAG IF NOT DOWN
16A9 3F 78         1770      CLR  SPEED1
16AB 3F 7C         1771      CLR  OUTSPD
16AD 81            1772      OVER1 RTS      ;RETURN IF NO BUTTONS DOWN
1773      ;
1774      ;      Set flag to increase the rate
1775      ;
1776      1776      BSET 0,CWF      ;SET FLAG FOR RATE INCREASE
1777      1777      JSR  SPEED      ;GO GET SPEED AND UPDATE
1778      1778      RTS
1779      ;
1780      ;      Set flag to decrease the rate
1781      ;
1782      1782      CCW1 BSET 0,CCWF      ;SET FLAG FOR RATE DECREASE
1783      1783      JSR  SPEED      ;GO SET SPEED AND UPDATE
1784      1784      RTS
1785      ;
1786      ;      Speed control for touch panel
1787      ;
1788      1788      SPEED INC  SPEED1
1789      1789      LDA  SPEED1

```

| | | | | |
|------------|------|----------|---------|----------------------------|
| 16BE A1 40 | 1790 | CMP | #\$40 | ;WAIT 1.98 SEC |
| 16C0 23 0A | 1791 | BLS | SPD1 | ;AT FIRST SPEED |
| 16C2 A1 90 | 1792 | CMP | #\$90 | ;WAIT 3.465 SEC |
| 16C4 23 15 | 1793 | BLS | SPD2 | |
| 16C6 A6 F0 | 1794 | LDA | #\$F0 | |
| 16C8 B7 78 | 1795 | STA | SPEED1 | |
| 16CA 20 21 | 1796 | BRA | SPD4 | |
| | 1797 | | | |
| | 1798 | | Speed 1 | |
| | 1799 | | | |
| 16CC 3C 7C | 1800 | SPD1 INC | OUTSPD | |
| 16CE B6 7C | 1801 | LDA | OUTSPD | |
| 16D0 A1 0F | 1802 | CMP | #\$0F | ;CHANGE EVERY .495 SEC |
| 16D2 27 19 | 1803 | BEQ | SPD4 | |
| 16D4 B6 78 | 1804 | LDA | SPEED1 | ;EXCEPT FIRST TIME THROUGH |
| 16D6 A1 01 | 1805 | CMP | #\$01 | |
| 16D8 27 13 | 1806 | BEQ | SPD4 | |
| 16DA 81 | 1807 | RTS | | |
| | 1808 | | | |
| | 1809 | | Speed 2 | |
| | 1810 | | | |
| 16DB 3C 7C | 1811 | SPD2 INC | OUTSPD | |

```

16DD B6 7C      1812      LDA      OUTSPD
16DF A1 06      1813      CMP      #$06      ;CHANGE EVERY .198 SEC
16E1 27 0A      1814      BEQ      SPD4
16E3 81         1815      RTS
1816      ;
1817      ;      Speed 3
1818      ;
16E4 3C 7C      1819      SPD3 INC      OUTSPD
16E6 B6 7C      1820      LDA      OUTSPD
16E8 A1 05      1821      CMP      #$05
16EA 27 01      1822      BEQ      SPD4
16EC 81         1823      RTS
1824      ;
16ED 3F 7C      1825      SPD4 CLR      OUTSPD
16EF 16 6D      1826      BSET 3,PSTAT1      ;MAKE START ON
16F1 18 6D      1827      BSET 4,PSTAT1      ;TURN DISPLAY ON
16F3 00 77 06   1828      BRSET 0,CCWF,DECL      ;IF CCW FLAG SET DEC DOSE OR RATE
1829      ;      ;ELSE INCREASE DOSE OR RATE
1830      ;      INCREASE
1831      ;
16F6 00 71 09   1832      BRSET 0,DFLAG,INCR1D      ;IF DOSE, INCREMENT DOSE
16F9 CC 131E    1833      JMP      INCRI      ;ELSE, INCREMENT RATE

```

SDOCID: <EP_____0327209A2_I_>

| | | | |
|--------------|------|-------------------|----------------------------------|
| 171B B7 55 | 1856 | STA DS1 | |
| 171D B6 56 | 1857 | LDA DS2 | |
| 171F A9 00 | 1858 | ADC #500 | ;BINARY ADD CARRY |
| 1721 8D | 1859 | DAA | ;CONVERT TO BCD |
| 1722 B7 56 | 1860 | STA DS2 | |
| 1724 CC 1745 | 1861 | JMP DSEXP | ; UPDATE DOSE |
| | 1862 | | |
| | 1863 | | |
| | 1864 | | |
| 1727 B6 57 | 1865 | DECRID LDA DOSECT | ;DELAY INC TO ALLOW DOSE DISPLAY |
| 1729 A1 1F | 1866 | CMP #31 | |
| 172B 23 1D | 1867 | BLS NODOSE | ;RETURN FROM INC |
| 172D A6 3D | 1868 | LDA #61 | ;RESET DOSE TIMER |
| 172F B7 57 | 1869 | STA DOSECT | ;EACH TIME BUTTON IS PUSHED |
| | 1870 | | |
| 1731 B6 55 | 1871 | LDA DS1 | ;IF DOSE IS 0 DONT DECREMENT |
| 1733 BA 56 | 1872 | ORA DS2 | |
| 1735 27 13 | 1873 | BEQ NODOSE | |
| | 1874 | | |
| | 1875 | | |
| 1737 B6 55 | 1876 | LDA DS1 | |
| 1739 AB 95 | 1877 | ADD #95 | ;10 COMPLIMENT BINARY SUBTRACT |

ISDOCID: <EP_____0327209A2_1_>


```

1900 ;
1901 ;
1902 ;      Update the digit numbers
1903 ;
1904 ;
1905      UPDATE      LDX  DSPDG1
1906                  LDA DIGITS,X
1907      STA  SAMPL4
1908      LDX  DSPDG2
1909      LDA  DIGITS,X
1910      STA  SAMPL3
1911      LDX  DSPDG3
1912      LDA  DIGITS,X
1913      STA  SAMPL2
1914      LDX  DSPDG4
1915      LDA  DIGITS,X
1916      STA  SAMPL1
1917      ;CHECK THIS CODE FOR HANDLING ZERO ST FLAG
1918 ;
1919      BCLR 6,ZEROST      ;INITIALIZE ZERO ST
1920      BCLR 7,ZEROST      ;AND
1921      CLR  COUNT3      ;COUNT3 FOR LED MUX ROUTINE

```

```

1922
1923 ;
1924 CLR HLDIM1 ;CLEAR HLDIM1 AND HLDIM2
1925 CLR HLDIM2 ;TO RESTART 2 1/2 MIN TIMER
1926 ;EACH TIME DISPLAY IS UPDATED
1927 ;
1928 ;
1929 ;
1930 ;
1931 CLR QH
1932 CLR PH
1933 LDA #10
1934 STA QL
1935 LDA DSPDG2
1936 STA PL
1937 JSR MULT16
1938 LDA DSPDG1
1939 ADD QL
1940 STA ALGO
1941 LDA #100
1942 STA QL
1943 LDA DSPDG3

```

CALCULATE MOTOR TIMES

```

177C 3F 67
177E 3F 69
1780 A6 0A
1782 B7 68
1784 B6 52
1786 B7 6A
1788 CD 196F
178B B6 51
178D BB 68
178F B7 46
1791 A6 64
1793 B7 68
1795 B6 53

```

| | | |
|--------------|------|------------|
| 1797 B7 6A | 1944 | STA PL |
| 1799 CD 196F | 1945 | JSR MULT16 |
| 179C B6 46 | 1946 | LDA ALGO |
| 179E B7 6A | 1947 | STA PL |
| 17A0 CD 1962 | 1948 | JSR ADD16 |
| 17A3 B6 67 | 1949 | LDA QH |
| 17A5 B7 69 | 1950 | STA PH |
| 17A7 B6 68 | 1951 | LDA QL |
| 17A9 B7 6A | 1952 | STA PL |
| 17AB A6 D6 | 1953 | LDA #\$D6 |
| 17AD B7 67 | 1954 | STA QH |
| 17AF A6 E8 | 1955 | LDA #\$E8 |
| 17B1 B7 68 | 1956 | STA QL |
| 17B3 CD 1993 | 1957 | JSR DIV16 |
| 17B6 B6 53 | 1958 | LDA DSPDG3 |
| 17B8 26 09 | 1959 | BNE MULT3 |
| 17BA B6 65 | 1960 | LDA XH |
| 17BC B7 48 | 1961 | STA TIMHI |
| 17BE B6 66 | 1962 | LDA XL |
| 17C0 B7 47 | 1963 | STA TIMLO |
| 17C2 81 | 1964 | RTS |
| | 1965 | ; |

| | | | | |
|--------------|------|---------|--------|----|
| 17C3 B6 65 | 1966 | MULT3 | IDA | XH |
| 17C5 B7 67 | 1967 | STA | QH | |
| 17C7 B6 66 | 1968 | IDA | XL | |
| 17C9 B7 68 | 1969 | STA | QL | |
| 17CB A6 03 | 1970 | IDA | #3 | |
| 17CD B7 6A | 1971 | STA | PL | |
| 17CF 3F 69 | 1972 | CIR | PH | |
| 17D1 CD 196F | 1973 | JSR | MULT16 | |
| 17D4 B6 68 | 1974 | IDA | QL | |
| 17D6 B7 47 | 1975 | STA | TIMIO | |
| 17D8 B6 67 | 1976 | IDA | QH | |
| 17DA B7 48 | 1977 | STA | TIMHI | |
| 17DC 81 | 1978 | RTS | | |
| | 1979 | ; | | |
| | 1980 | ; | | |
| | 1981 | \$EJECT | | |

```

1982 ;
1983 ;
1984 ; AC OR DC CHECK
1985 ; -----
1986 ;
1987 ACDC BRCIR 6,PORTC,AC ;IF CLR ITS AC
1988 BCLR 1,ACON ;IF NOT CLEAR THE FLAG
1989 RTS
1990 AC BSET 1,ACON ;SET FLAG ITS AC
1991 ; CLR T25 ;RESET TIMER FOR AUTO OFF
1992 ; CLR T20
1993 RTS
1994 ;
1995 $EJECT

```

```

1996 ;
1997 ;
1998 ; BATTERY CHECK
1999 ; -----
2000 ;
2001 17E6 02 6E 37 BATCK BRSET 1, ACON, NOLOB ; IF ON AC DONT DO
2002 17E9 00 80 45 BRSET 0, DSPTST, NLB ; DONT DO IF TEST MODE
2003 17EC 00 6E 1D BRSET 0, LOWBAT, LBIM ; IF LOW ALREADY TIME 15 MIN.
2004 17EF 03 6D 07 BRCIR 1, PSTAT1, LOBCK ; IF MOTOR OFF CHECK
2005 17F2 B6 5B LDA REG1 ; WAIT UNTILL MOTOR RUNS .4 SEC
2006 17F4 A1 1E CMP #30
2007 17F6 22 01 BHI LOBCK
2008 17F8 81 RTS
2009 17F9 08 02 35 LOBCK BRSET 4, PORTC, NLB ; IF NO SIGNAL THEN LOW BAT
2010 17FC 05 6D 1A BRCIR 2, PSTAT1, KILTIM ; IF PUMP IS OFF, DISABLE PUMP
2011 17FF 10 6E BSET 0, LOWBAT ; SET LOW BAT FLAG
2012 1801 1E 6E BSET 7, ERRON ; SET ERROR
2013 1803 1A 6E BSET 5, ALRM ; SET ALARM
2014 1805 17 6D BCIR 3, PSTAT1 ; GET OUT OF START MODE TO ENABLE LOW BAT
2015 1807 11 71 BCIR 0, DFLAG ; CLEAR DOSE MODE FLAG
2016 1809 1D 6E BCIR 6, VOIM ; CLEAR VOLUME MODE FLAG
2017 180B 81 RTS

```

```

2018 ;
2019 ; INCREMENT 15 MIN TIMER
2020 ;
2021 LBIM INC BATTM1 ;INCREMENT BATTERY TIMER: BATTM1,BATTM2
2022 BNE CONCNT
2023 INC BATTM2
2024 CONCNT LDA BATTM2
2025 CMP #$07 ;IF BAT TIMER HAS COUNTED 15 MIN
2026 BEQ KILTIM ;KILL POWER
2027 RTS
2028 KILTIM SEI ;SET INT MASK TO PREVENT TURN-ON
2029 LDA #$0F ;DISABLE PUMP
2030 STA PORTB
2031 BRA KILTIM
2032 ;
2033 NOLOB BRCIR 0,LOWBAT,NLB ;DONT CLEAR ANY THING IF NO LOW BAT
2034 ;THIS MAINTAINS OTHER ERRORS IF PRESENT
2035 BCIR 0,LOWBAT ;CLR LOW BAT FLAG
2036 BCIR 7,ERRON ;CLR ERROR FLAG
2037 BCIR 5,ALRM ;CLR ALARM
2038 BSET 4,PSTAT1 ;MAKE SURE DISPLAY ON
2039 JSR UPDATE ;UPDATE DISPLAY

```

SDOCID: <EP_____0327209A2_1_>


```

2049 ;
2050 ;      DEAD BATTERY CHECK
2051 ;
2052 ;      This routine checks the dead battery signal (portc(7)).
2053 ;      If this signal is active (portc(7)=0), then the routine
2054 ;      kills power to the processor.
2055 ;-----
2056 ;
2057 ;      DEATCK      BRSET 1,ACON,NODECK      ;DON'T DO IF AC POWERED
2058 ;      BRCLR 1,PSTAT1,DEADCK ;IF MOTOR OFF CHECK
2059 ;      LDA REG1      ;WAIT UNTILL MOTOR RUNS .4 SEC
2060 ;      CMP #30
2061 ;      BHI DEADCK
2062 ;      RTS
2063 ;
2064 ;      DEADCK BRSET 7,PORTC,NODECK ;IF NO DEAD BAT SIGNAL, RETURN
2065 ;
2066 ;      DBKIL SEI      ;SET INT MASK TO PREVENT TURN-ON
2067 ;      LDA #0F      ;DISABLE PUMP
2068 ;      STA PORTB
2069 ;      BRA DBKIL      ;LOOP UNTIL POWER DISSIPATES
2070 ;

```

| | 184D 81 | 2071 | NODECK | RTS |
|--|---------|------|---------|-----|
| | | 2072 | ; | |
| | | 2073 | \$EJECT | |

```
2074 ;  
2075 ; ALARM  
2076 ;-----  
2077 ; Conditions  
2078 ;  
2079 ALARM BRCLR 2,PSTAT1,NOAIM ;IF OFF DONT ALARM  
2080 BRSET 6,VOLM,NOAIM ;IF IN VOLUME DONT DO  
2081 BRSET 5,ALRM,ALMRT ;IF ALARM FLAG SET THEN ALARM  
2082 ;  
2083 ;  
2084 ALCLR BCCLR 7,PORTB ;SET ALARM HIGH  
2085 NOAIM RTS  
2086 ALMRT BRSET 4,PSTAT1,ALCLR  
2087 BSET 7,PORTB  
2088 RTS  
2089 ;  
2090 $EJECT
```

SDOCID: <EP_____0327209A2_1_>

7593

```

187E 04 6E 50      2113      BRSET 2,NDROP,FLO1      ;IF FLOW ERR THEN DO
1881 00 70 45      2114      BRSET 0,DOSE1,DOSE1      ;IF DOSE ERROR THEN DO
1884 06 6E 46      2115      BRSET 3,HOLDER,HLD1      ;IF HOLD ERROR THEN DO
1887 20 59          2116      BRA ERR1
1889 81             2117      NOBLIN      RTS
2118      ;
2119      ;          Hold mode blink
2120      ;
188A 18 6D          2121      SETBL BSET 4,PSTAT1      ;BLANK DISPLAY
188C 0F 6E FA      2122      BRCLR 7,ERRON,NOBLIN
188F 00 6E 54      2123      BRSET 0,LOWBAT,BAT1      ;GO DO BAT IF LO BAT
1892 08 6E 5E      2124      BRSET 4,NOSET,SET2      ;GO DO SET IF NO SET
1895 04 6E 4A      2125      BRSET 2,NDROP,ERR1      ;GO DO ERR IF FLOW ERR
1898 00 70 26      2126      BRSET 0,DOSE1,OUT2      ;GO DO DEL IF DOSE ERR
189B 06 6E 44      2127      BRSET 3,HOLDER,ERR1      ;GO DO ERR IF HOLD ERR
189E 20 3E          2128      BRA SYS1
18A0 81             2129      RTS
2130      ;
2131      ;          Blink the DP during the run mode
2132      ;
18A1 BE 63          2133      DOTBL LDH COUNT2      ;GET DP #
18A3 E6 4D          2134      LDA SAMPLE1,X      ;GET NEW DIGIT

```

```

18A5 AA 80      2135      ORA   #$80      ;CLEAR DP
18A7 E7 4D      2136      STA   SAMPL1,X  ;STORE IT FOR DISPLAY
18A9 5C         2137      INCX          ;BUMP POINTER TO NEXT DIGIT
18AA A3 04      2138      CPX   #$04      ;DONE ALL THREE
18AC 27 09      2139      BEQ   NEWX      ;RESET IF SO
18AE E6 4D      2140      DOTON LDA SAMPL1,X ;GET NEXT DIGIT
18B0 A4 7F      2141      AND   #$7F      ;TURN DP ON
18B2 E7 4D      2142      STA   SAMPL1,X  ;STORE IT FOR DISPLAY
18B4 BF 63      2143      STX   COUNT2    ;SAVE COUNT
18B6 81         2144      RTS
                ;
18B7 0F 6D 03   2145      NEWX BRC1R 7,PSTAT1,DOT24 ;IF 224, SCROLL 3 DIGITS
18BA 5F         2146      C1RX          ;RESET COUNTER
18BB 20 F1      2147      BRA   DOTON    ;GO TURN ON
                ;
18BD AE 01      2148      DOT24 LDX   #$01    ;IF THIRD DIGIT
18BF 20 ED      2149      BRA   DOTON    ;ELSE CONTINUE WITH BLINK ROUTINE
                ;
                ;      Error messages
                ;
18C1 A6 8C      2150      OUT2  LDA   #$8C
18C3 B7 4D      2151      STA   SAMPL1
                ;

```

```

2157 ;
18C5 AE 16 LDX #$16 ;DISPLAY ?
2158 BRA OUTCHR
18C7 20 2E
2159 ;
2160 DOSE1 LDX #$19 ;DISPLAY DOSE
18C9 AE 19 BRA OUTCHR
18CB 20 2A
2161 ;
2162 ;
2163 HLD1 LDX #$12 ;POINT TO ERR MESSAGE
2164 BRA OUTCHR ;GO LOAD UP DIGITS
18CD AE 12
2165 ;
2166 FLO1 LDX #$06 ;POINT TO FLO MESSAGE
2167 BRA OUTCHR ;GO LOAD UP DIGITS
18D1 AE 06
2168 ;
2169 ;
2170 ;
18D5 AE 00 L01 LDX #$00
18D7 CD 18F7 JSR OUTCHR
18DA 01 6D C4 BRC1R 0, PSTAT1, DOTBL ;IF IN RUN MODE, UPDATE DP SCROLL
18DD 81 RTS
2171 ;
2172 SYS1 LDX #$1C ;DISPLAY SYS
18DE AE 1C
2173 ;
2174 ;
2175 ;
2176 ;
2177 ;
2178 ;

```

```
2179      ;
2180      ERR1  LDX  #$09
2181      18E4 20 11  BRA  OUTCHR
2182      ;
2183      18E6 AE 03  BAT1  LDX  #$03
2184      18E8 CD 18F7 JSR  OUTCHR
2185      18EB 01 6D B3  BRC1R 0,PSTAT1,DOTBL  ;IF IN RUN MODE, UPDATE DP SCROLL
2186      18EE 81      RTS
2187      ;
2188      18EF AE 0C      NO1  LDX  #$0C
2189      18F1 20 04      BRA  OUTCHR
2190      ;
2191      18F3 AE 0F      SET2  LDX  #$0F
2192      18F5 20 00      BRA  OUTCHR
2193      ;
2194      18F7 D6 1000  OUTCHR  LDA  TABLE5,X
2195      18FA B7 4E      STA  SAMPLE2
2196      18FC 5C      INCX
2197      18FD D6 1000  LDA  TABLE5,X
2198      1900 B7 4F      STA  SAMPLE3
2199      1902 5C      INCX
2200      1903 D6 1000  LDA  TABLE5,X
```


| | | | |
|------------|------|---------|--------|
| 1906 B7 50 | 2201 | STA | SAMPL4 |
| 1908 81 | 2202 | RTS | |
| | 2203 | ; | |
| | 2204 | \$EJECT | |

7593

```

2205 ;
2206 ; HOLD ERROR
2207 ; -----
2208 ; Conditions
2209 ;
2210 HILDER BRCIR 2,PSTAT1,NOHE ;DONT DO IF OFF
2211 BRCIR 0,PSTAT1,NOHE ;DONT DO IF NOT IN HOLD
2212 ;
2213 ; Test for hold error
2214 ;
2215 INC HLDTIM1 ;INCREMENT HOLD COUNTER
2216 BNE CNCNT
2217 INC HLDTIM2
2218 CNCNT LDA HLDTIM2 ;CHECK IF 2 1/2 MIN ELAPSED
2219 CMP #501
2220 BHI HILDER1 ;IF HLDTIM2 > 1, HOLDER
2221 BEQ CONCK ;IF HLDTIM2= 1, CHECK HLDTIM1
2222 RTS ;ELSE RETURN
2223 ;
2224 CONCK LDA HLDTIM1
2225 CMP #527 ;IF HLDTIM1 >= 27, HOLDER
2226 BHS HILDER1

```

```
1924 81      2227      RTS
                2228      ;
1925 17 6D    2229      HLDER1      BCLR 3,PSTAT1      ;CLEAR START FLAG
1927 16 6E    2230      BSET 3,HOLDER      ;SET HOLD ERROR FLAG
1929 1A 6E    2231      BSET 5,ALRM      ;ALARM ON
192B 1E 6E    2232      BSET 7,ERRON      ;ERROR FLAG ON
192D 81    2233      NOHE RTS
                2234      ;
                2235      $EJECT
```

```

2236 ;
2237 ; ONE DAY WAIT FOR AUTO OFF
2238 ;-----
2239 ; Conditions
2240 ;
2241 TIM24 BRSET 2,PSTAT1,KEEPON ;IF UNIT ON DONT DO
2242 ;
2243 ; 24HR TURN OFF
2244 ;
2245 INC T25 ;GET COUNTER
2246 LDA T25
2247 CMP #SEF ;(EF + 1) = 240 x 6 = 24 HOURS
2248 BLO KEEPON
2249 ;
2250 KILLP SEI ;SET INT MASK TO PREVENT TURN-ON
2251 LDA #0F ;DISABLE PUMP
2252 STA PORTB
2253 BRA KILLP ;LOOP UNTIL POWER DISSAPTES
2254 ;
2255 KEEPON RTS
2256 ;
2257 ;

```

192E 04 6D 0F

1931 3C 45

1933 B6 45

1935 A1 EF

1937 25 07

1939 9B

193A A6 0F

193C B7 01

193E 20 F9

1940 81

```
2258 ; CLEAR VOLUME SUBROUTINE
2259 ; This routine is called by tim24 and ckclr routines.
2260 ; -----
2261 ;
2262 CLVOL CLR DECML1
2263 CLR VS1
2264 CLR VS2
2265 RTS
2266 ;
2267 ; CLEAR TIMER COUNTERS
2268 ; -----
2269 ;
2270 TIMRST CLR TS
2271 CLR T5
2272 CLR T10
2273 CLR T15
2274 RTS
2275 ; COMPARE DOSE TO VOLUME
2276 ; -----
2277 ;
2278 ; IF VOLUME DELIVERED >= PROGRAMMED DOSE, SET X REGISTOR TO FF
2279 ;
```

| | | | | |
|------------|------|--------|-----------|----------------------------|
| 1951 5F | 2280 | CPVTOD | CIRX | |
| 1952 B6 5A | 2281 | | LDA VS2 | |
| 1954 B1 56 | 2282 | | CMP DS2 | |
| 1956 25 09 | 2283 | | BLO VLTD | ;IF VOLUME < DOSE, CLEAR C |
| 1958 22 06 | 2284 | | BHI VGTD | ;IF VOLUME > DOSE, SET C |
| | 2285 | | | ;ELSE, CHECK VS1 |
| | 2286 | | LDA VS1 | |
| 195A B6 59 | 2287 | | CMP DS1 | |
| 195C B1 55 | 2288 | | BLO VLTD | ;IF VOLUME < DOSE, CLEAR C |
| 195E 25 01 | 2289 | | | |
| | 2290 | | VGTD COMX | ;SET X TO FF |
| 1960 53 | 2291 | | VLTD RTS | |
| 1961 81 | 2292 | | | |
| | 2293 | | | |
| | 2294 | | \$EJECT | |

```

2295 ;*****
2296 ;
2297 ;
2298 ;          MATH UTILITIES FOLLOW
2299 ;
2300 ;
2301 ;*****
2302 ;
2303 ;
2304 ;*****
2305 ;
2306 ;  PROGRAM ADDS 2, 16-BIT UNSIGNED BINARY NUMBERS, PRODUCING A 17-BIT
2307 ;  RESULT.
2308 ;          ENTER WITH:  2, UNSIGNED 16-BIT OPERANDS TO BE ADDED IN
2309 ;                      PH , PL and QH , QL.
2310 ;          EXIT WITH:   17-BIT RESULT IN:  CARRY , QH, QL
2311 ;                      (QH, QL DESTROYED) .
2312 ;
2313 ;*****
2314 ;
2315 ;
2316 ADD16 LDA  QL      ;ADD 16 BYTES.

```

```

1964 BB 6A      2317      ADD      PL
1966 B7 68      2318      STA      QL
1968 B6 67      2319      LDA      QH      ;ADD MS BYTES.
196A B9 69      2320      ADC      PH
196C B7 67      2321      STA      QH      ;16-BIT RESULT IN QH, QL. OVERFLOW IN CARRY.
                2322      ;
196E 81         2323      RTS
                2324      ;
                2325      ;
                2326      ;
                2327      ;
                2328      ;*****
                2329      ;      *
                2330      ;      PROGRAM MULTIPLIES 2, 16 BIT UNSIGNED BINARY OPERANDS, CREATING A 32--*
                2331      ;      BIT UNSIGNED RESULT. (NO OVERFLOW IS POSSIBLE).      *
                2332      ;      ENTER WITH:  OPERANDS TO BE MULTIPLIED IN:      *
                2333      ;                      QH , QL      *
                2334      ;                      and PH , PL      *
                2335      ;      EXIT WITH:  32-BIT RESULT IN:      XH , XL , QH , QL      *
                2336      ;                      *
                2337      ;*****
                2338      ;

```



```

2339      ;
2340      MULT16 LDX      #16      ;LOOP COUNTER.
2341      CLR      XH      ;CLEAR UPPER 16 BITS OF 32-BIT ACCUM.
2342      CLR      XL
2343      ROR      QH      ;CHECK BIT 0 OF QL.
2344      ROR      QL
2345      NXT      BCC      ROTAT  ;IF 0, DON'T ADD, JUST SHIFT.
2346      LDA      XL      ;OTHERWISE, ADD IN THE CONTENTS OF PH , PL TO
2347      ADD      PL      ;XH , XL.
2348      STA      XL
2349      LDA      XH
2350      ADC      PH
2351      STA      XH
2352      ;
2353      ROTAT      ROR      XH      ;SHIFT THE 32-BIT ACCUM. 1 BIT RIGHT.
2354      ROR      XL
2355      ROR      QH
2356      ROR      QL
2357      ;
2358      DECX
2359      BNE      NXT
2360      ;

```

196F AE 10
1971 3F 65
1973 3F 66
1975 36 67
1977 36 68
1979 24 0C
197B B6 66
197D BB 6A
197F B7 66
1981 B6 65
1983 B9 69
1985 B7 65

1987 36 65
1989 36 66
198B 36 67
198D 36 68

198F 5A
1990 26 E7

1992 81
2361 RTS ;OTHERWISE, RETURN WITH RESULT IN XH,XL,QH,QL.
2362 ;
2363 ;
2364 ;*****
2365 ; *
2366 ; PROGRAM PERFORMS THE DIVISION OF 2, 16 BIT UNSIGNED OPERANDS, PRODUC-
2367 ; ING A 16 BIT UNSIGNED RESULT: *
2368 ; *
2369 ; (QH , QL/ PH , PL) -----> XH , XL *
2370 ; *
2371 ; ENTER WITH: 16 BIT DIVISOR IN PH , PL *
2372 ; 16 BIT DIVIDEND IN QH , QL *
2373 ; *
2374 ; EXIT WITH: QUOTIENT TRUNCATED TO 16 BITS *
2375 ; IN XH , XL *
2376 ; REGISTERS AFFECTED: X, A, COUNT1, TEMPX, TEMPY *
2377 ; (QH, QL, PH, PL DESTROYED) *
2378 ; *
2379 ;*****
2380 ;
2381 ;
2382 DIV16 LDA #1
1993 A6 01

| | | | | |
|------------|------|-------------|--------|---|
| 1995 3D 69 | 2383 | TST | PH | |
| 1997 2B 0B | 2384 | BMI | DIV153 | ;IF DIVISOR IS LEFT-JUSTIFIED. |
| 1999 4C | 2385 | DIV151 INCA | | ;OTHERWISE, KEEP SHIFTING DIVISOR LEFT |
| 199A 38 6A | 2386 | ASL | PL | ;UNTIL THE MSB IN PH = 1, OR UNTIL |
| 199C 39 69 | 2387 | ROL | PH | ;16 SHIFTS HAVE BEEN DONE. |
| 199E 2B 04 | 2388 | BMI | DIV153 | |
| 19A0 A1 11 | 2389 | CMP | #17 | |
| 19A2 26 F5 | 2390 | BNE | DIV151 | |
| 19A4 B7 62 | 2391 | STA | COUNT1 | ;COUNT1 = # SHIFTS REQUIRED +1. |
| 19A6 B6 67 | 2392 | LDA | QH | ;MOVE THE DIVIDEND INTO A, X. |
| 19A8 BE 68 | 2393 | LDX | QL | |
| 19AA 3F 67 | 2394 | CIR | QH | ;MAKE WAY FOR THE QUOTIENT. |
| 19AC 3F 68 | 2395 | CIR | QL | |
| 19AE BF 6C | 2396 | STX | TEMPX | ;STORAGE FOR THE DIVIDEND AFTER SUBTRACTING |
| 19B0 B7 6B | 2397 | STA | TEMPA | ;OUT DIVISOR. |
| 19B2 9F | 2398 | TXA | | |
| 19B3 B0 6A | 2399 | SUB | PL | ;TRY SUBTRACTING THE DIVISOR. |
| 19B5 B7 6C | 2400 | STA | TEMPX | |
| 19B7 B6 6B | 2401 | LDA | TEMPA | ;SAVE THE REMAINDER IN TEMPA, TEMPX. |
| 19B9 B2 69 | 2402 | SBC | PH | |
| 19BB B7 6B | 2403 | STA | TEMPA | |
| 19BD BE 6C | 2404 | LDX | TEMPX | |

| | | | | | |
|------------|------|------------|--------|--------|---|
| 19BF 24 10 | 2405 | | BCC | DIV165 | ;IF CARRY=0, THEN DIVISOR WAS SMALLER THAN |
| | 2406 | | | | ;DIVIDEND. GO SET THE CURRENT QUOTIENT BIT. |
| | 2407 | | | | ;OTHERWISE, ADD THE DIVISOR BACK IN, |
| 19C1 9F | 2408 | TXA | | | |
| 19C2 BB 6A | 2409 | ADD | PL | | |
| 19C4 B7 6C | 2410 | STA | TEMPX | | |
| 19C6 B6 6B | 2411 | LDA | TEMPA | | |
| 19C8 B9 69 | 2412 | ADC | PH | | |
| 19CA B7 6B | 2413 | STA | TEMPA | | ;AND SAVE IT. |
| 19CC BE 6C | 2414 | LDX | TEMPX | | |
| 19CE 98 | 2415 | CLC | | | ;THE QUOTIENT BIT WILL BE 0. |
| 19CF 20 01 | 2416 | BRA | DIV167 | | |
| 19D1 99 | 2417 | DIV165 SEC | | | ;THE QUOTIENT BIT WILL BE 1. |
| 19D2 39 68 | 2418 | DIV167 ROL | QL | | ;ROTATE THE QUOTIENT LEFT 1 BIT, |
| 19D4 39 67 | 2419 | ROL | QH | | ;SHIFTING THE MOST RECENT QUOTIENT BIT |
| 19D6 34 69 | 2420 | LSR | PH | | ;INTO THE LSB. |
| 19D8 36 6A | 2421 | ROR | PL | | |
| 19DA 3A 62 | 2422 | DEC | COUNT1 | | ;KEEP GOING UNTIL THE COUNTER=0. |
| 19DC 26 D0 | 2423 | BNE | DIV163 | | |
| 19DE B6 67 | 2424 | LDA | QH | | ;WHEN DONE, MOVE THE RESULT INTO XH, XL. |
| 19E0 B7 65 | 2425 | STA | XH | | |
| 19E2 B6 68 | 2426 | LDA | QL | | |

```

19E4 B7 66      2427      STA      XL
                  ;
19E6 81      2428
                2429      RTS      ;RETURN.
                2430
                2431      ;*****
                2432      ;
                2433      ;          BCDEXP
                2434      ;
                2435      ;  PROGRAM CONVERTS 2 BYTE BCD NUMBER POINTED TO BY THE X REGISTER TO *
                2436      ;  A 4 BYTE DECIMAL NUMBER. THE FOUR BYTE NUMBER IS CONVERTED ON THE FLY *
                2437      ;  TO THE CORRECT LED DISPLAY SEGMENT CODE WHICH IS SENT TO THE DISPLAY *
                2438      ;  BY THE LED MUX ROUTINE.
                2439      ;
                2440      ;          ENTER WITH: ADDRESS OF LOWER BCD BYTE IN X
                2441      ;
                2442      ;          EXIT WITH: LED DISPLAY CODE FOR 4 DECIMAL DIGITS
                2443      ;          IN SAMPLE1, SAMPLE2, SAMPLE3, SAMPLE4
                2444      ;          *
                2445      ;          REGISTERS AFFECTED: X, A, SAMPLE1 - SAMPLE4, TEMPX
                2446      ;
                2447      ;*****
                2448

```

```

19E7 BF 6C      2449      BCDEXP      STX  TEMPX      ;STORE POINTER FOR LATER USE
19E9 F6         2450      LDA  ,X      ;GET LOWER BYTE OF BCD DIGIT
19EA A4 0F      2451      AND  #$0F    ;MASK UPPER NIBBLE
19EC 97         2452      TAX
19ED D6 101F    2453      LDA  DIGITS,X ;CONVERT DECIMAL TO LED CODE
19F0 B7 50      2454      STA  SAMPL4
19F2 BE 6C      2455      ;
19F4 F6         2456      LDX  TEMPX    ;RESTORE FIRST BCD BYTE
19F5 44         2457      LDA  ,X
19F6 44         2458      LSR4
19F7 44         2459      LSR4
19F8 44         2460      LSR4
19F9 97         2461      LSR4
19FA D6 101F    2462      TAX
19FD B7 4F      2463      LDA  DIGITS,X ;CONVERT DECIMAL TO LED CODE
19FF 3C 6C      2464      STA  SAMPL3
1A01 BE 6C      2465      ;
1A03 F6         2466      INC  TEMPX    ;SET POINTER TO UPPER BCD DIGIT
1A04 A4 0F      2467      LDX  TEMPX    ;GET SECOND BCD DIGIT
1A06 97         2468      LDA  ,X
1A08 97         2469      AND  #$0F    ;MASK UPPER NIBBLE
1A0A 97         2470      TAX

```

```

1A07 D6 101F      2471      LDA  DIGITS,X      ;CONVERT DECIMAL TO LED CODE
1A0A B7 4E         2472      STA  SAMPL2
                    2473      ;
1A0C BE 6C         2474      LDX  TEMPX          ;RESTORE SECOND BCD DIGIT
1A0E F6           2475      LDA  ,X
1A0F 44           2476      LSR  A
1A10 44           2477      LSR  A
1A11 44           2478      LSR  A
1A12 44           2479      LSR  A          ;PUSH UPPER NIBBLE TO LOWER NIBBLE
1A13 97           2480      TAX
1A14 D6 101F      2481      LDA  DIGITS,X      ;CONVERT DECIMAL TO LED CODE
1A17 B7 4D         2482      STA  SAMPL1
                    2483      ;
                    2484      ;      UPDATE DISPLAY FLAGS USED BY LED MUX
                    2485      ;
1A19 1D 74        2486      BCIR 6,ZEROST      ;CLEAR MSD ZERO FLAG
1A1B 1F 74        2487      BCIR 7,ZEROST      ;CLEAR MSD-1 ZERO FLAG
1A1D 3F 64        2488      CTR  COUNT3          ;CLEAR MUX COUNTER
1A1F 81           2489      RTS
                    2490      ;
                    2491      SEJECT

```

| | | | |
|------|------------|------------|---------|
| 2492 | | | |
| 2493 | =1C00 | ORG | PMP1ST |
| 2494 | | | |
| 2495 | 1C00 1E 01 | ELOOP BSET | 7,PORTB |
| 2496 | 1C02 20 FC | BRA | ELOOP |
| 2497 | | | ; |
| 2498 | | | END |

FIG. 1

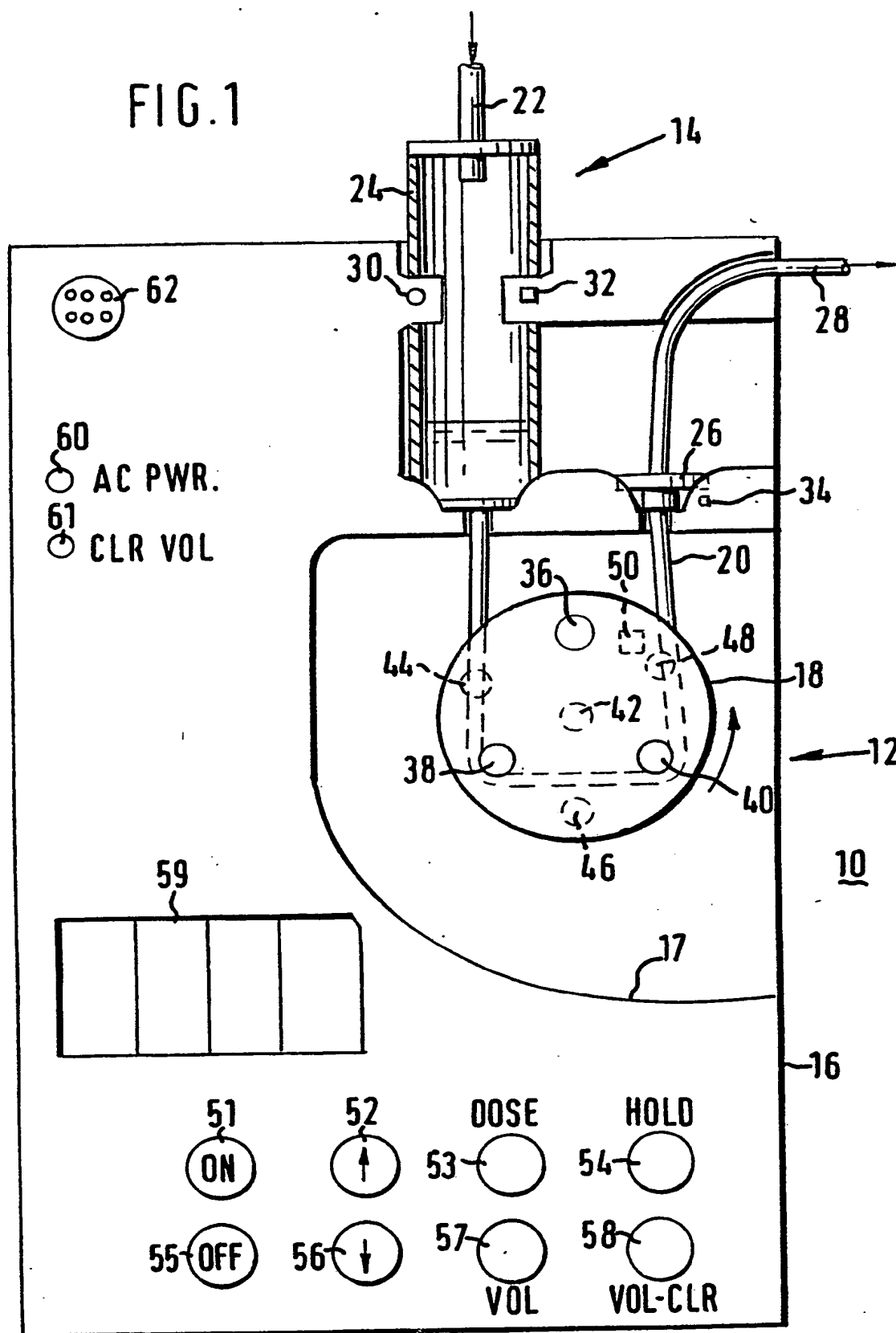


FIG. 2a

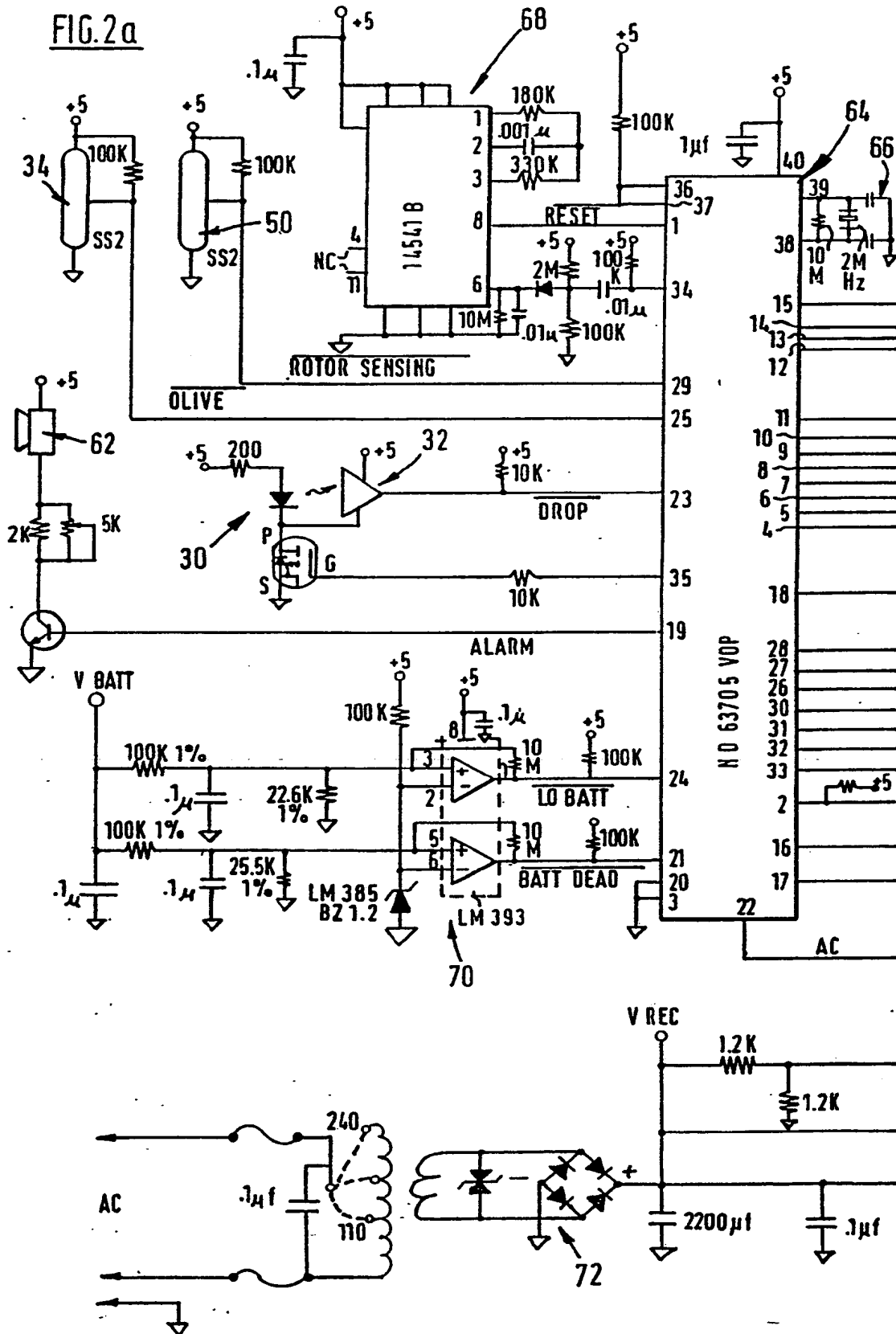


FIG. 2b

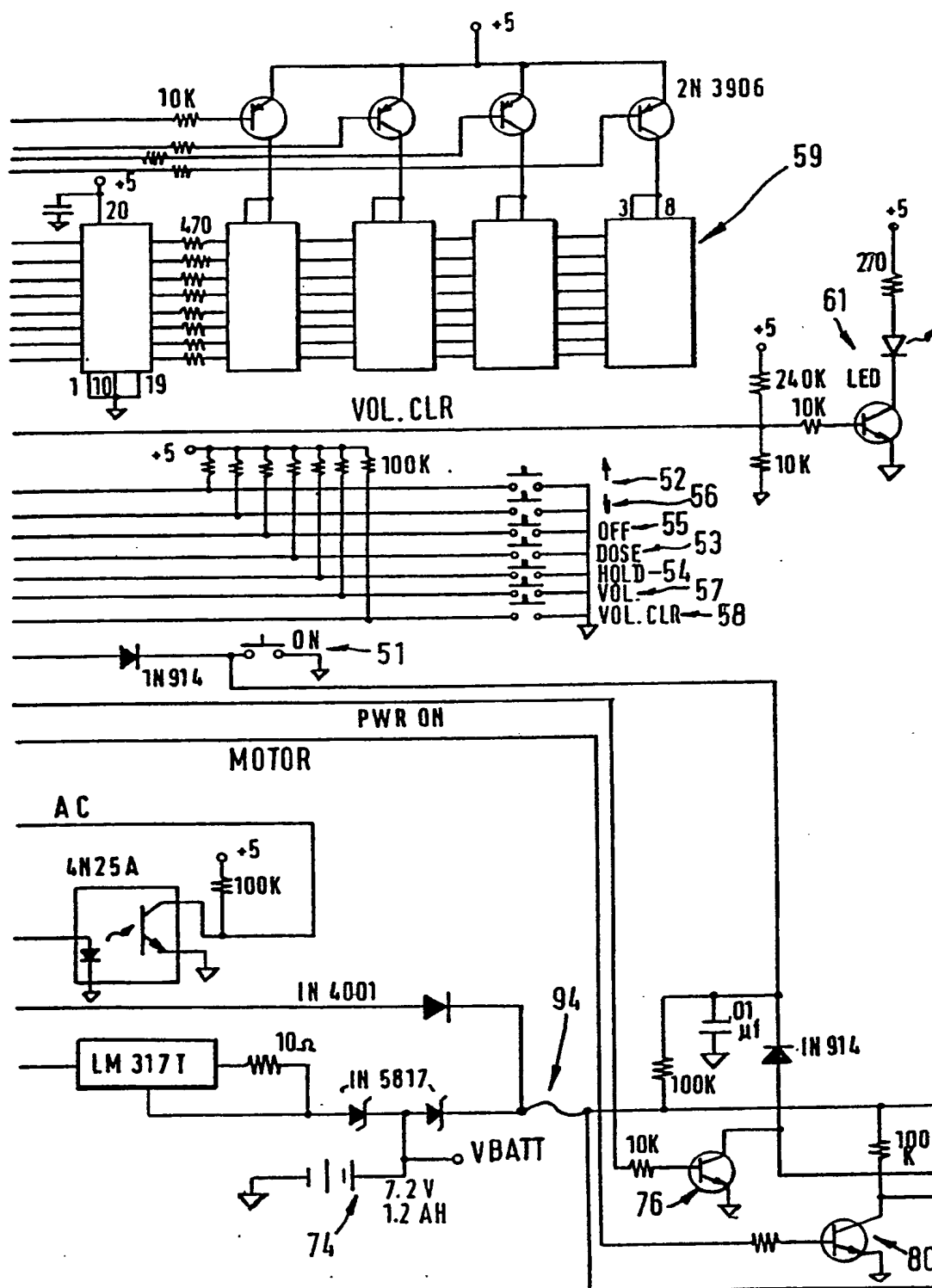
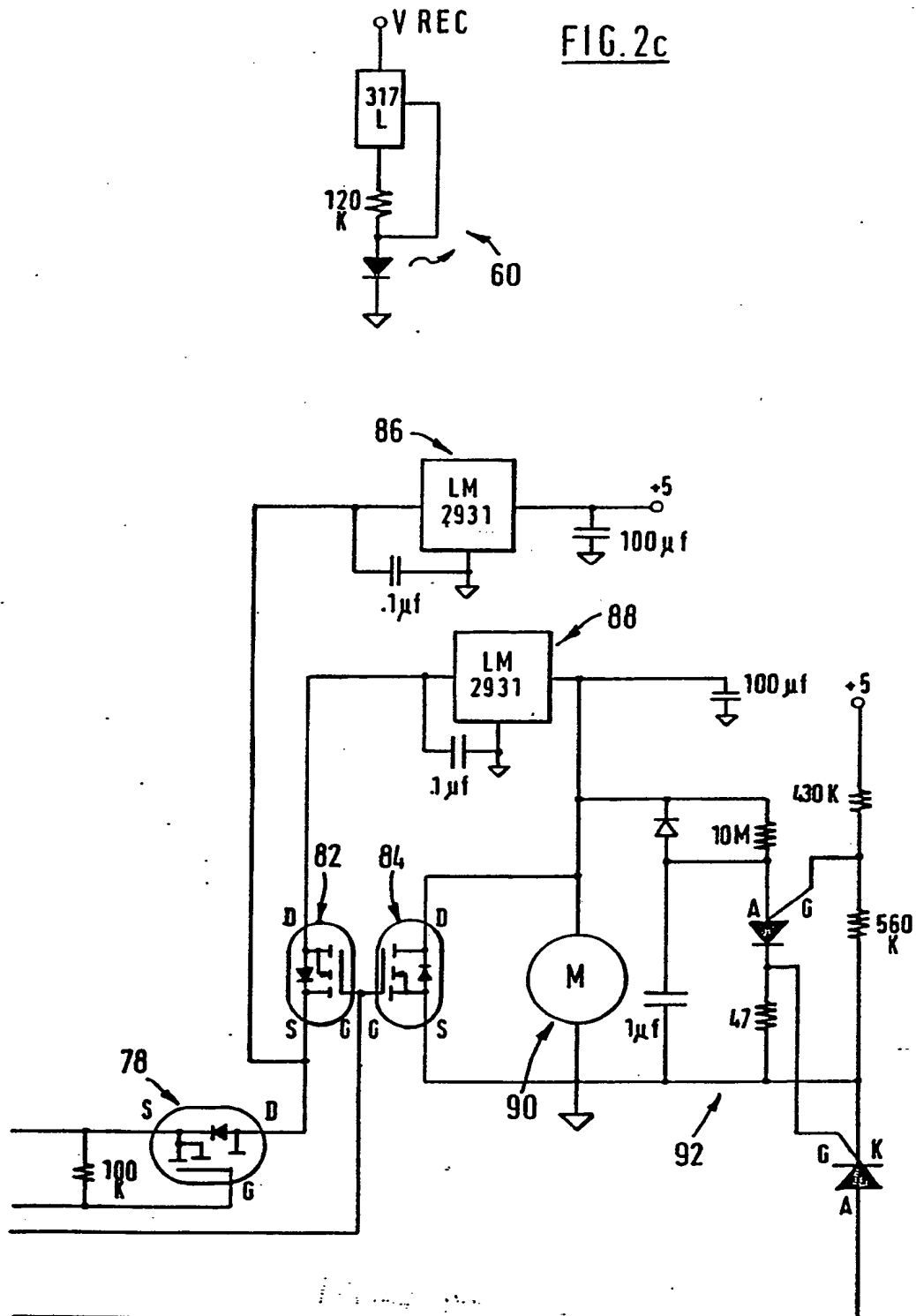


FIG. 2c



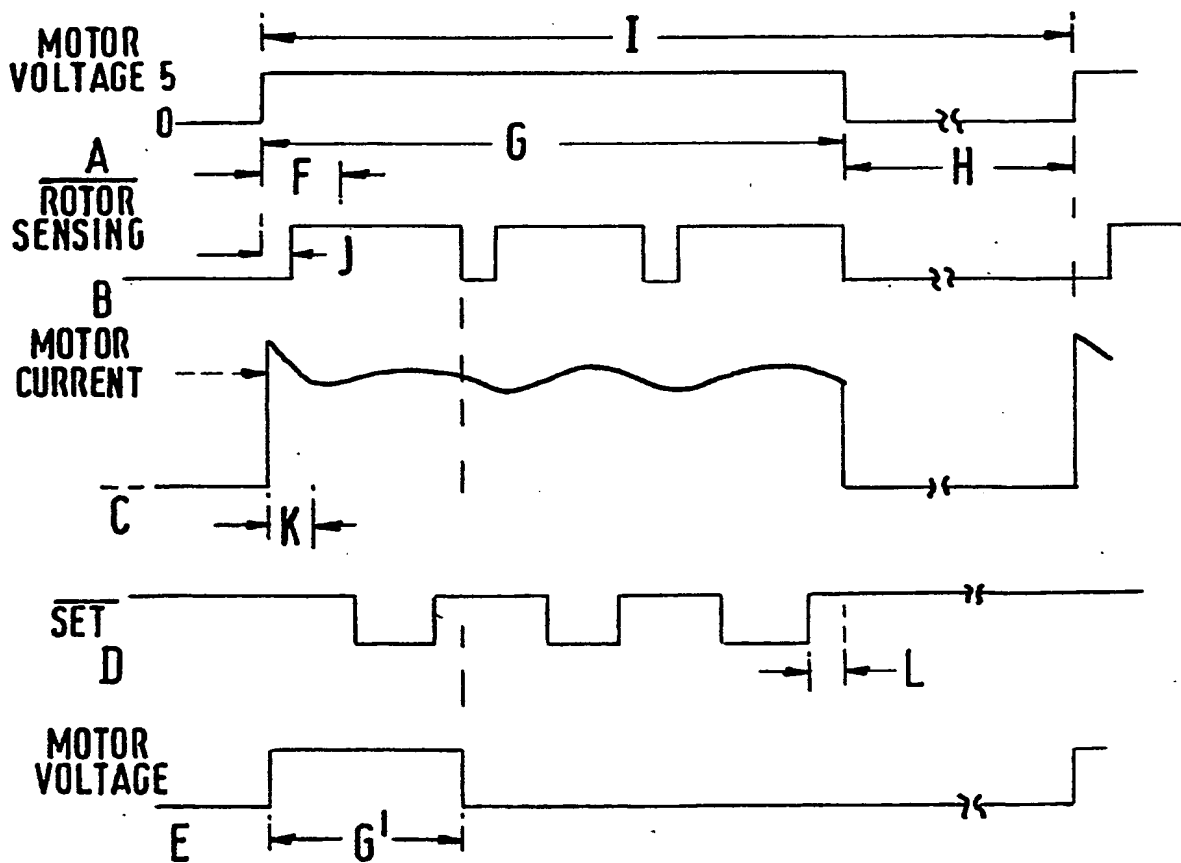
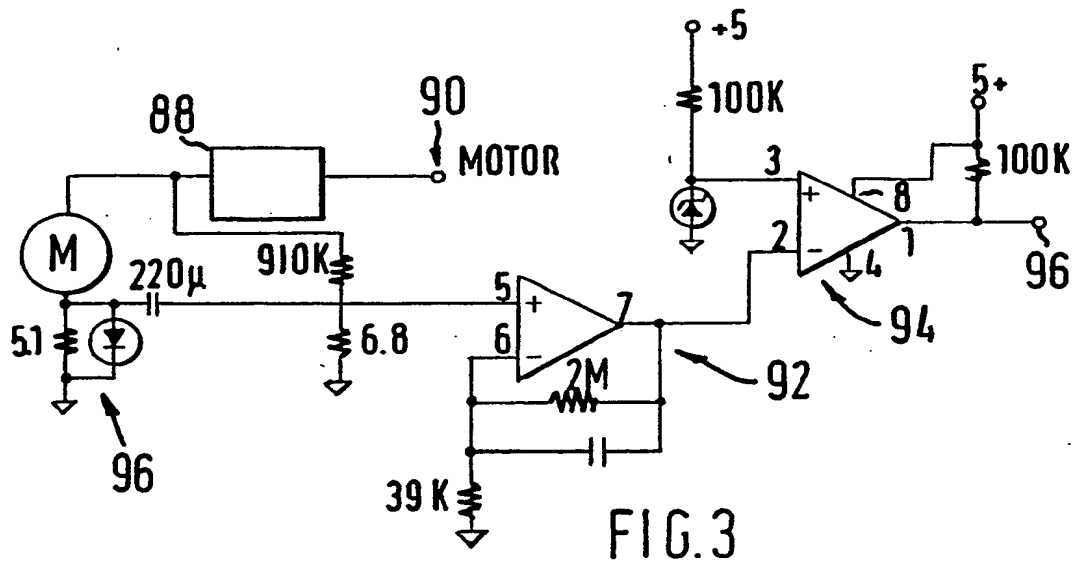


FIG 4

⑫

EUROPEAN PATENT APPLICATION

⑪ Application number: 89300296.4

⑤① Int. Cl.⁴: **A 61 M 5/14**
F 04 B 51/00

⑫ Date of filing: 13.01.89

③① Priority: 15.01.88 US 144795

④③ Date of publication of application:
09.08.89 Bulletin 89/32

⑥④ Designated Contracting States:
AT BE CH DE ES FR GB GR IT LI LU NL SE

⑧⑧ Date of deferred publication of search report:
15.11.89 Bulletin 89/46

⑦① Applicant: **SHERWOOD MEDICAL COMPANY**
1831 Olive Street
St. Louis, MO 63103 (US)

⑦② Inventor: **Jackson, Edward**
18 Gail Drive
Northford Connecticut 06472 (US)


Pasqualucci, Joseph
12 Spruce Brook Road
Seymour Connecticut 06483 (US)

⑦④ Representative: **Brown, Keith John Symons et al**
c/o John Wyeth & Brother Limited Patent and Trademark
Department Huntercombe Lane South
Taplow Maidenhead Berkshire SL6 0PH. (GB)

⑥④ **Motor unit for a fluid pump and method of operation.**

⑤⑦ An enteral nutrition pump system (10) operates in a cyclical manner with a period between cycles being selected in accordance with the desired fluid delivery rate. Each pump cycle may correspond to a single rotation of the rotor (18) or a fractional rotation of the rotor. Rotor rotation may alternatively be sensed by utilization of magnetic sensors (50) or by monitoring of the AC component of current supplied to a DC motor driving the rotor.

FIG.1

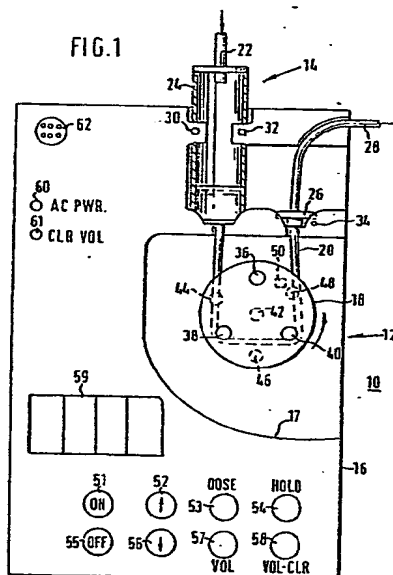

Bundesdruckerei Berlin



EP 89 30 0296

| DOCUMENTS CONSIDERED TO BE RELEVANT | | | |
|---|---|--|--|
| Category | Citation of document with indication, where appropriate, of relevant passages | Relevant to claim | CLASSIFICATION OF THE APPLICATION (Int. Cl. 4) |
| X | DE-A-2 651 962 (SIEMENS AG) * Page 10, line 22 - page 11, line 5, line 25 - page 12, line 15; figure 3 * | 1-4,6,8-11 | A 61 M 5/14 F 04 B 51/00 |
| Y | -- | 12 | |
| X | EP-A-0 127 346 (PERITRONIC) * Page 8, line 10 - page 10, line 22 * | 1-3,8,9 | |
| Y | -- | 12 | |
| X | US-A-4 498 843 (SCHNEIDER et al.) * Abstract; column 9, line 44 - column 11, line 34 * | 1-3,8,9 | |
| X | WO-A-86 01 413 (SCHWEIZER) * Claims 1-5; abstract * | 1-3 | TECHNICAL FIELDS SEARCHED (Int. Cl. 4) A 61 M H 02 P F 04 B |
| X | GB-A-2 011 652 (NIKKISO CO. LTD) * Page 3, lines 36-75 * | 1,2 | |
| A | EP-A-0 090 152 (SGS-ATES) * Abstract; figure 4 * | 5,7,13 | |
| A | US-A-3 610 779 (HUBBY) * Column 4, line 62 - column 5, line 71; figures 4,5 * | 5,7,13,15 | |
| The present search report has been drawn up for all claims | | | |
| Place of search THE HAGUE | | Date of completion of the search 02-08-1989 | Examiner CLARKSON |
| CATEGORY OF CITED DOCUMENTS | | | |
| X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document | | T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons & : member of the same patent family, corresponding document | |

EPO Form 1503 03/82



CLAIMS INCURRING FEES

The present European patent application comprised at the time of filing more than ten claims.

- ☐ All claims fees have been paid within the prescribed time limit. The present European search report has been drawn up for all claims.
- ☐ Only part of the claims fees have been paid within the prescribed time limit. The present European search report has been drawn up for the first ten claims and for those claims for which claims fees have been paid.
namely claims:
- ☐ No claims fees have been paid within the prescribed time limit. The present European search report has been drawn up for the first ten claims.

X LACK OF UNITY OF INVENTION

The Search Division considers that the present European patent application does not comply with the requirement of unity of invention and relates to several inventions or groups of inventions.

namely:

1. Claims 1-4,6,8-12: Pump control involving magnetic sensing
2. Claims 1-3,5,7,13-15: Pump control involving current sensing of DC motor

- ☒ All further search fees have been paid within the fixed time limit. The present European search report has been drawn up for all claims.
- ☐ Only part of the further search fees have been paid within the fixed time limit. The present European search report has been drawn up for those parts of the European patent application which relate to the inventions in respect of which search fees have been paid.
namely claims:
- ☐ None of the further search fees has been paid within the fixed time limit. The present European search report has been drawn up for those parts of the European patent application which relate to the invention first mentioned in the claims.
namely claims:

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☒ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.